



US 20250277278A1

(19) **United States**(12) **Patent Application Publication**
DEAN et al.(10) **Pub. No.: US 2025/0277278 A1**(43) **Pub. Date: Sep. 4, 2025**(54) **QUALITATIVE AND QUANTITATIVE
DETERMINATION OF SINGLE VIRUS
HAPLOTYPES IN COMPLEX SAMPLES**(30) **Foreign Application Priority Data**

Apr. 9, 2020 (EP) 20169080.7

(71) Applicant: **Takeda Vaccines, Inc.**, Cambridge, MA
(US)**Publication Classification**(72) Inventors: **Hansi DEAN**, Weston, MA (US);
Lovkesh KARWAL, Newton, MA
(US); **Jill A. LIVENGOOD**,
Cambridge, MA (US); **Subash DAS**,
Cambridge, MA (US); **Callum BELL**,
Santa Fe, NM (US); **Johnny SENA**,
Santa Fe, NM (US); **Nicholas
DEVITT**, Santa Fe, NM (US)(51) **Int. Cl.**
C12Q 1/70 (2006.01)
C12Q 1/686 (2018.01)
C12Q 1/6869 (2018.01)
(52) **U.S. Cl.**
CPC **C12Q 1/70** (2013.01); **C12Q 1/686**
(2013.01); **C12Q 1/6869** (2013.01); **C12Q**
2600/172 (2013.01)(73) Assignee: **Takeda Vaccines, Inc.**, Cambridge, MA
(US)(57) **ABSTRACT**


The present invention provides methods for detecting the presence of at least one revertant virus in a sample containing an attenuated virus by single-molecule real-time (SMRT) sequencing. It further provides methods for the quantitative analysis of a mixture of virus haplotypes in a virus sample by SMRT sequencing. It also provides kits useful in performing said methods.

Specification includes a Sequence Listing.(21) Appl. No.: **17/911,152**(22) PCT Filed: **Apr. 9, 2021**(86) PCT No.: **PCT/US2021/026530**

§ 371 (c)(1),

(2) Date: **Sep. 13, 2022**

PacBio Amplicon Sequencing Workflow

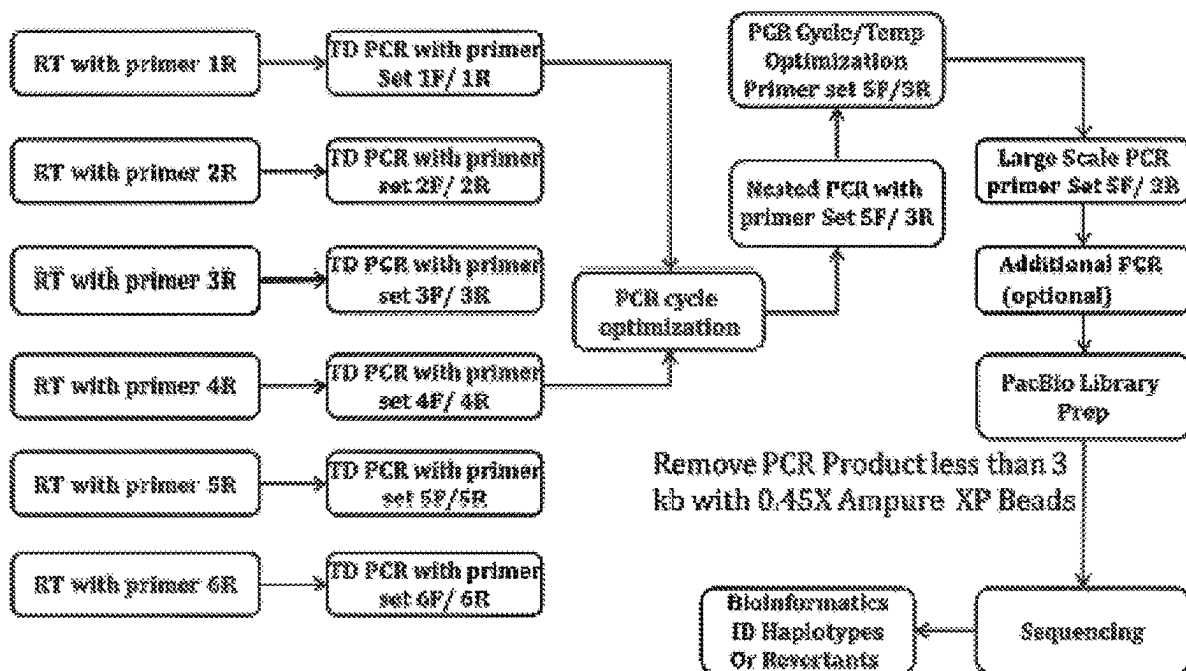

PacBio Amplicon Sequencing Workflow

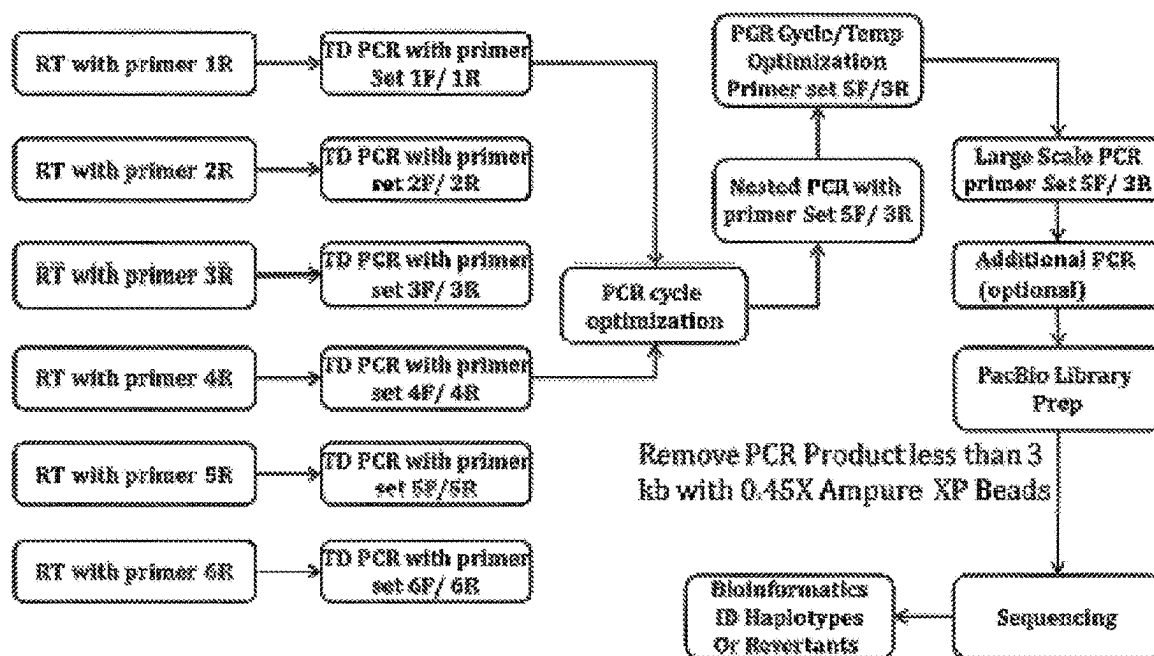


FIG. 1

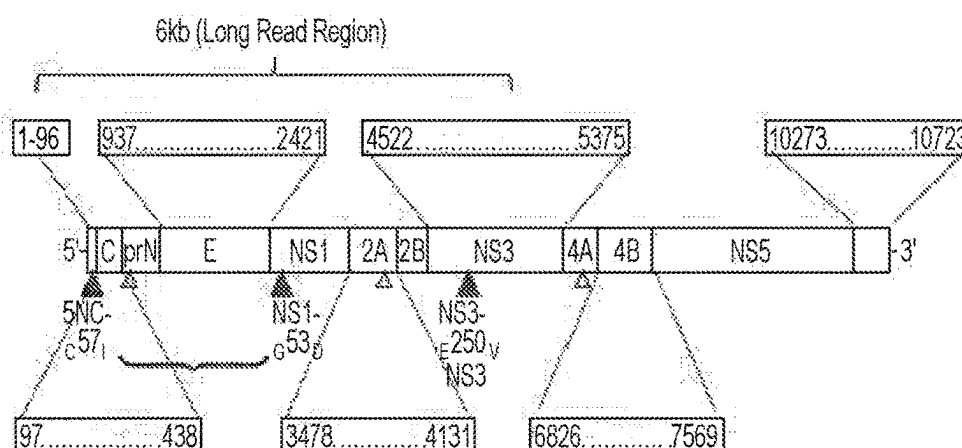


FIG. 2

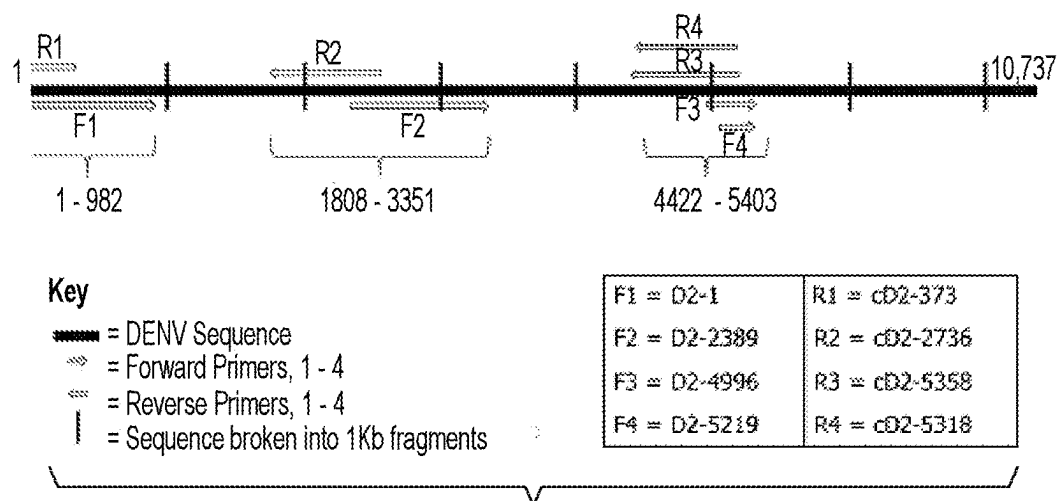


FIG. 3

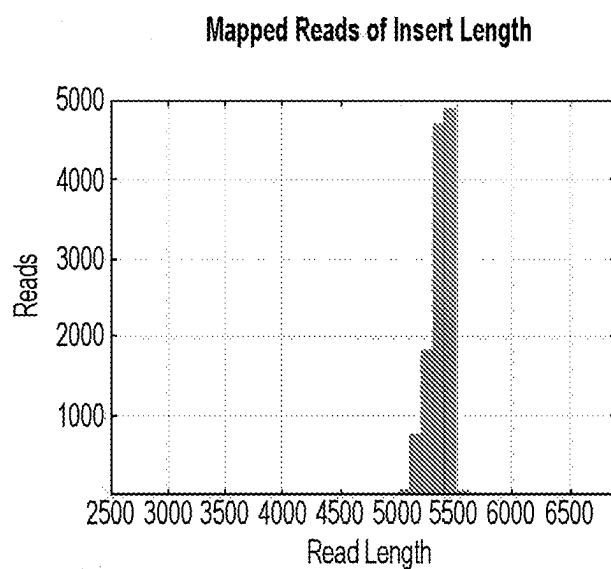


FIG. 4

QUALITATIVE AND QUANTITATIVE DETERMINATION OF SINGLE VIRUS HAPLOTYPES IN COMPLEX SAMPLES

CROSS REFERENCE TO RELATED APPLICATIONS

[0001] This application is a 35 U.S.C. § 371 National Stage filing of International Application No. PCT/US2021/026530, filed Apr. 9, 2021, which claims the benefit of priority to European Patent Application No. 20169080.7, filed Apr. 9, 2020. The entire contents of the aforementioned applications are incorporated herein by reference in their entireties.

SEQUENCE LISTING

[0002] This application contains a Sequence Listing which is submitted herewith in electronically readable ASCII format via EFS-Web, and is hereby incorporated by reference in its entirety. The electronic Sequence Listing was created on Nov. 23, 2022, and was filed on Mar. 24, 2023, and is named “2022-12-23_Corrected Sequence Listing-TXT” and is 20 KB in size.

TECHNICAL FIELD

[0003] The present invention relates to methods for detecting the presence of at least one revertant virus in a sample containing an attenuated virus. It also relates to methods for the quantitative analysis of a mixture of virus haplotypes in a virus sample. It further relates to kits useful in performing said methods.

BACKGROUND OF THE INVENTION

[0004] Dengue virus is an arbovirus in the Flaviviridae family. It is an enveloped virus with an isometric capsid and a linear positive sense RNA genome. Dengue viruses are endemic in most regions of the tropics and are transmitted by infectious mosquitoes, typically *Aedes aegypti*, causing symptomatic infections or asymptomatic seroconversion. Symptomatic dengue infection is a systemic and dynamic disease. It has a wide clinical spectrum that includes both severe and non-severe clinical manifestations. There are four different dengue serotypes (DENV 1-4) which are genetically related to yellow fever and similar tick-borne encephalitis viruses.

[0005] One tetravalent dengue vaccine (TDV) presently tested in clinical trials consists of all four dengue serotypes, TDV1, TDV2, TDV3 and TDV4. TDV2 is a live attenuated recombinant virus that originated from PDK-53 that became non-pathogenic after continuous passaging of DEN-2 virus in primary dog kidney (PDK) cells. The mutations necessary and sufficient for the attenuated phenotype of TDV2 virus have been genetically identified (Butrapet et al., J. Virol., vol. 74, no. 7 (2000), 3011-3019). These mutations are C-to-T at genomic nucleotide position 57 in the 5' noncoding region, and Gly-to-Asp and Glu-to-Val mutations at amino acid positions 53 and 250 in nonstructural proteins (NS) 1 and 3, respectively. This corresponds to a C to T mutation at genomic nucleotide position 57, a G- to A mutation at genomic nucleotide position 2579 and an A to T mutation at genomic nucleotide position 5270. These three genetic determinants of attenuation reside outside of the structural gene regions of the TDV2 viral genome. TDV1, TDV3 and TDV4 are chimeric variants having the backbone of Dengue

virus 2 PDK53 and the prM and E proteins of Dengue virus 1, Dengue virus 3 and Dengue virus 4, respectively, which maintain the attenuation phenotype of TDV-2.

[0006] TDV is comprised of four viruses when administered into humans. Host immune response is triggered due to vaccine virus replication. Amplification of virus may generate mutations leading to reversion to wild type. In a scenario when there are febrile cases in the first 30 days post vaccination, it is important to confirm that these illnesses are not caused by vaccine viruses that may have been reverted to wild type and became pathogenic, but may be due to an infection by other closely related Flaviviruses or any other microbial agents. In such a complicated scenario it is important to decipher and identify the genome variants present in the suspected sample. This requires sequencing methods that can accurately capture the genetic diversity of RNA virus genomes.

[0007] WO 2017/179017 A1 describes that for sensitive detection of the presence of revertants in samples containing attenuated viruses Tagman-based mismatch amplification mutation assay (Taq-MAMA) was used which is not based on sequencing.

[0008] Huang et al., PLOS Negl. Trop. Diseases 7 (2013), e2243, 1-11 relates to the genetic and phenotypic characterization of manufacturing seeds of a tetravalent Dengue vaccine (DENVax). For sequencing of the virus Sanger sequencing on an automated sequencer was used. For testing of the reversion at position 5'NC-57 locus of the vaccine a Tagman-based mismatch amplification mutation assay (Taq-MAMA) was used. This represents an improved quantitative SNP assay. According to the authors the detection sensitivity could be lowered to detect 0.01-0.07% reversion. The authors further report that the reversion sensitivity limit detectable by consensus genome sequence analysis is 10 to 30%. The person skilled in the art would therefore have considered that genome sequence analysis is not suitable for sensitive detection of revertants in samples containing attenuated viruses.

[0009] In addition, possible revertants have been tested in a murine neurovirulence model. The prior art therefore suffers from the drawbacks that no sequence information of entire single virus haplotypes such as revertants can be provided. Due to the lack of identification of entire single virus haplotypes the prior art cannot quantify the amount of entire single virus haplotypes in complex mixtures such as patient samples or live attenuated vaccines.

[0010] Ardui et al., Nucl. Acids Res. 46 (2018), 2159-2168 describes that SMRT sequencing has been performed in hepatitis B virus, hepatitis C virus and HIV-1. None of those studies involved the detection of revertants in samples containing attenuated viruses.

[0011] Dilernia et al., Nucl. Acids Res. 43 (2015), e129, pages 1 to 13 referred to by Ardui et al. describes the application of a complex algorithm used for error correction of the SMRT sequencing method. As can be seen from Table 1, genomic variants having an expected frequency of less than 1.56% could not be detected at all. For higher frequencies, the observed value at least to some extent differed significantly from the expected value. The investigated genomes differed in at least five nucleotides from each other.

[0012] In view thereof, there is a need for the provision of improved methods for the detection of revertants in complex samples containing live attenuated viruses, in particular Dengue virus. Further, there is a need for the quantitation of

single viral haplotypes in complex samples, in particular Dengue virus containing samples.

SUMMARY OF THE INVENTION

[0013] The technical problems underlying the invention are solved by the provision of the subject-matter as defined in the claims.

[0014] According to a first aspect, a method for detecting the presence of at least one revertant virus in a sample containing an attenuated virus comprising the steps:

[0015] a) purifying virus nucleic acids from said sample;

[0016] b) optionally, preparing double-stranded DNA templates from said virus nucleic acids;

[0017] c) amplifying said viral nucleic acids of a) or said double-stranded DNA templates of b) by polymerase chain reaction (PCR) such that PCR products having a length of at least 2 kb and comprising at least one attenuation locus of said attenuated virus are generated;

[0018] d) sequencing said PCR products by single-molecule real-time (SMRT) sequencing; and

[0019] e) comparing at least one nucleotide at said at least one attenuation locus in said PCR products with the corresponding nucleotide in the wild-type virus, wherein the presence of the wild-type nucleotide at said at least one attenuation locus is indicative of the presence of a revertant virus in said sample.

[0020] The inventors have found that SMRT sequencing (PacBio sequencing) is useful for the qualitative and quantitative determination of revertants within a sample of live attenuated virus vaccine. SMRT sequencing was known to exhibit a high sequencing error rate of 13% e.g. in comparison to a low sequencing error rate of 0.001% of the conventional Sanger method. Due to the use of high number or Read counts and additionally high number of UMI counts the PacBio inherent high sequence error is significantly lowered. The successful use of SMRT sequencing for the determination of single or multiple revertants in a sequence having a length of more than 5 kb, wherein the revertants occur in low quantities in a live attenuated vaccine sample may be considered as a surprising advantageous effect. The use of SMRT sequencing provides the advantage of a high read length which is useful for covering the entire Dengue virus genome. A further advantage of SMRT sequencing of live attenuated vaccine samples over known approaches such as three separate PCR sequencing reactions is that the entire sequence of a single haplotype can be determined in comparison to an “average” haplotype when conventional PCR approaches are used to determine the sequence at the attenuation loci. An average haplotype means that the sequence is not determined by sequencing of a single nucleic acid molecule, but by sequencing of different nucleic acid molecules which cover different parts of the complete nucleic acid molecule. Hence, the nucleic acid sequence determined at the second or third attenuation locus may not be derived from the same nucleic acid molecule used for determining the sequence at the first attenuation locus. While Sanger sequencing is useful for determining the presence of a revertant in a sample comprising different nucleic acid molecules, it cannot determine the whole nucleic acid sequence of the revertant in one sequencing reaction.

[0021] SMRT sequencing was known to exhibit lower throughput compared to second generation sequencing technologies. This was considered an obstacle in quantitative analysis, such as gene/isoform abundance estimation. To overcome this obstacle third generation sequencing technologies such as SMRT may be evaluated as hybrid sequencing, i.e. a combination of second generation and third generation sequencing technology data. It was therefore surprising that the use of SMRT sequencing alone was sufficient to provide a quantitative analysis of complex mixtures of viruses, such as Dengue viruses. The present inventors in preferred embodiments have provided a highly sensitive method for the detection of revertants in samples containing attenuated viruses. The quantity of the revertant virus which may be detected may be as low as 0.2% based on the total amount of virus in the sample. This is significantly better than prior art approaches for HIV virus.

[0022] In addition, in preferred embodiments of the present invention at least two different haplotypes present in the sample can be distinguished from each other which may only differ in a single nucleotide from each other. This is significantly better than prior art approaches for HIV virus which could only distinguish haplotypes being at least five nucleotides different from each other.

[0023] According to a second aspect, a method for the quantitative analysis of a mixture of virus haplotypes in a virus sample is provided comprising the steps:

[0024] a) purifying virus nucleic acids from said sample;

[0025] b) optionally, preparing double-stranded DNA templates from said virus nucleic acids by using a primer comprising a unique molecular identifier (UMI) sequence such that said double-stranded DNA templates comprise the UMI sequence;

[0026] c) amplifying said double-stranded DNA templates by polymerase chain reaction (PCR) such that PCR products having a length of at least 2 kb and comprising a unique molecular identifier (UMI) sequence are generated;

[0027] d) sequencing said PCR products by single-molecule real-time (SMRT) sequencing; and

[0028] e) determining the relative amounts of said PCR products and thereby the virus haplotypes contained in the sample.

[0029] Surprisingly, the inventors have found that a UMI workflow using the above method identified 8 to 11 haplotypes per revertant, i.e. is more sensitive compared to a BLASR or Quiver workflow which only identified 2 to 5 haplotypes per revertant on the basis of the same sample. Samples of the Dengue virus reference monorevertants for reversion at positions NS1-2579 (P1), NS3-5270 (P3) and 5'NC-57 (P5 and P51), respectively, were obtained from Center of Disease Control and Prevention (CDC). In said reference samples the UMI workflow identified significantly more individual haplotypes compared to the BLASR or the Quiver workflow.

[0030] According to a third aspect the use of the method according to the first aspect in the quality control of vaccines containing live, attenuated virus is provided.

[0031] According to a fourth aspect the use of the method according to the second aspect in the diagnosis of a virus infection in a patient sample is provided.

[0032] According to a fifth aspect is provided a kit for detecting the presence of at least one revertant virus in a

sample containing an attenuated dengue virus comprising a first primer pair comprising the nucleotide sequence of SEQ ID NO: 9 (5F) and a second primer comprising the nucleotide sequence of SEQ ID NO: 6 (3R).

[0033] According to a sixth aspect is provided a kit for the quantitative analysis of a mixture of virus haplotypes in a dengue virus sample comprising a first primer comprising the nucleotide sequence of SEQ ID NO:21 (Primer IIA) and a second primer comprising the nucleotide sequence of SEQ ID NO:9 (5F).

BRIEF DESCRIPTION OF THE DRAWINGS

[0034] FIG. 1 is a chart showing an example of the SMRT (Pacbio) sequencing work flow. PCR means Polymerase Chain Reaction; RT means Reverse Transcription PCR; TD PCR means Touchdown PCR; ID haplotypes means identification of haplotypes; ID revertants means identification of revertants.

[0035] FIG. 2 shows the TDV2 genome.

[0036] FIG. 3 shows the positions of primers on TDV2 genome used for Sanger sequencing.

[0037] FIG. 4 shows the mapped reads of Insert Length.

DETAILED DESCRIPTION OF THE INVENTION

[0038] Where the term “comprise” or “comprising” is used in the present description and claims, it does not exclude other elements or steps. For the purpose of the present invention, the term “consisting of” is considered to be an optional embodiment of the term “comprising”. If hereinafter a group is defined to comprise at least a certain number of embodiments, this is also to be understood to disclose a group which optionally consists only of these embodiments.

[0039] Where an indefinite or a definite article is used when referring to a singular noun e.g. “a” or “an”, “the”, this includes a plural form of that noun unless specifically stated. Vice versa, when the plural form of a noun is used it refers also to the singular form.

[0040] Furthermore, the terms first, second, third or (a), (b), (c) and the like in the description and in the claims are used for distinguishing between similar elements and not necessarily for describing a sequential or chronological order. It is to be understood that the terms so used are interchangeable under appropriate circumstances and that the embodiments of the invention described herein are capable of operation in other sequences than described or illustrated herein.

[0041] In the context of the present invention any numerical value indicated is typically associated with an interval of accuracy that the person skilled in the art will understand to still ensure the technical effect of the feature in question. As used herein, the deviation from the indicated numerical value is in the range of $\pm 10\%$, and preferably of $\pm 5\%$. The aforementioned deviation from the indicated numerical interval of $\pm 10\%$, and preferably of $\pm 5\%$ is also indicated by the terms “about” and “approximately” used herein with respect to a numerical value.

[0042] According to the first aspect of the present invention is provided a method for detecting the presence of at least one revertant virus in a sample containing an attenuated virus comprising the steps:

[0043] a) purifying virus nucleic acids from said sample;

[0044] b) optionally, preparing double-stranded DNA templates from said virus nucleic acids;

[0045] c) amplifying said viral nucleic acids of a) or said double-stranded DNA templates of b) by polymerase chain reaction (PCR) such that PCR products having a length of at least 2 kb and comprising at least one attenuation locus of said attenuated virus are generated;

[0046] d) sequencing said PCR products by single-molecule real-time (SMRT) sequencing; and

[0047] e) comparing at least one nucleotide at said at least one attenuation locus in said PCR products with the corresponding nucleotide in the wild-type virus sequence, wherein the presence of the wild-type nucleotide at said at least one attenuation locus is indicative of the presence of a revertant virus in said sample.

[0048] “Virus” herein means any virus including double-stranded and single-stranded DNA viruses, and double and single-stranded RNA viruses.

[0049] “An attenuated virus” herein means that the virus has a reduced replication capacity and/or a reduced infectivity in host cells compared to wild-type virus. The replication capacity and/or infectivity may be determined in vitro in suitable cell systems or in vivo in suitable animal models. Attenuation may be achieved by serial passaging of the virus in a foreign host such as in tissue culture, embryonated eggs or live animals. Alternatively, attenuation may be performed by chemical agents.

[0050] Live attenuated viruses are important for use in viral vaccines, since the live attenuated virus generates a stronger immune response compared to an inactivated virus. The live attenuated virus may be an RNA virus or a DNA virus. Suitable live attenuated RNA viruses may be selected from dengue virus, poliovirus, rubella virus, measles virus, yellow fever virus, mumps virus and influenza virus. Preferably, the live attenuated RNA virus is dengue virus. Suitable live attenuated DNA viruses include varicella zoster virus, vaccinia virus, smallpox virus, Herpes simplex virus and chikungunya virus. The live attenuated virus may differ from the wild-type virus in one or more mutations in the nucleic acid sequence of the virus. For example, for one tetravalent dengue vaccine (TDV) it is known that attenuation of neurovirulence in newborn mice is determined by mutations 5'NCR-57 C->U, NS1-53 Gly->Asp (nt-2579 G->A) and NS3-250 Glu->Val (nt-5270 A->U) (Butrapet et al., J. Virol. 74 (2000), 3011-3019). The above nucleotide sequence numbering relates to the nucleotide sequence of the tetravalent Dengue virus serotype 2 (TDV2) set forth in SEQ ID NO:23.

[0051] “A revertant virus” herein means a virus derived from an attenuated virus having a reversion of a nucleotide at one or more attenuation loci to the wild-type sequence. The detection of revertant virus is of particular importance for the quality control of vaccines on the basis of live attenuated viruses. In the case of TDV, a revertant will have nucleotide C at position 57, nucleotide G at position 2579 and/or nucleotide A on position 5270.

[0052] “Sample” herein means any sample from vaccinated individuals or patients or a sample to be tested at the stage of quality control of manufactured vaccines. The sample may be from whole blood or serum, most preferred

the samples are from serum. The sample may also be vaccine composition or the stock solution prior to administration. This may be important for quality control and safety reasons.

[0053] Preferably a vaccine composition to be analysed according to the present invention comprises a dengue antigen of each of serotypes 1 to 4 which are each independently selected from the group consisting of: (a) a live attenuated dengue virus and (b) a live attenuated chimeric virus.

[0054] Preferably a vaccine composition to be analysed according to the present invention comprises a dengue antigen of each of serotypes 1 to 4, wherein said dengue antigens of serotypes 1, 3 and 4 are each a live attenuated chimeric dengue virus and said dengue antigen of serotype 2 is selected from the group consisting of a live attenuated dengue virus and a live attenuated chimeric dengue virus. For example, a vaccine composition to be analysed of the present invention may comprise a dengue antigen of each of serotypes 1 to 4, wherein said dengue antigens of serotypes 1, 3 and 4 are each a live attenuated chimeric dengue/dengue virus and said dengue antigen of serotype 2 is a live attenuated dengue virus. For example, a vaccine composition to be analysed according to the present invention may be the tetravalent mixture of dengue antigens of each of serotypes 1 to 4 (referred to as DENVax) which is disclosed in Huang et al., *PLoS Negl Trop Dis* 7(5): e2243 (2013). Alternatively, a vaccine composition to be analysed according to the present invention may comprise a dengue antigen of each of serotypes 1 to 4, wherein said dengue antigens of serotypes 1, 3 and 4 are each a live attenuated chimeric YF/dengue virus and said dengue antigen of serotype 2 is a live attenuated dengue virus.

[0055] Preferably a vaccine composition to be analysed according to the present invention comprises a dengue antigen of each of serotypes 1 to 4, wherein said dengue antigens of serotypes 1, 3 and 4 are each independently selected from the group consisting of: (a) a live attenuated dengue virus and (b) a live attenuated chimeric dengue virus and said dengue antigen of serotype 2 is a live attenuated chimeric dengue virus. For example, a vaccine composition to be analysed according to the present invention may comprise a dengue antigen of each of serotypes 1 to 4, wherein said dengue antigens of serotypes 1, 3 and 4 are each a live attenuated dengue virus and said dengue antigen of serotype 2 is a live attenuated chimeric dengue/dengue virus. For example, a vaccine composition according to the present invention may be any of the tetravalent mixtures of dengue antigens of each of serotypes 1 to 4 (referred to as TV001, TV002, TV003 and TV004') which are disclosed in Durbin et al., *Journal of Infectious Diseases* (2013), 207, 957-965. Preferably, a vaccine composition to be analysed according to this embodiment of the invention is TV003.

[0056] Preferably a vaccine composition to be analysed according to the present invention comprises a dengue antigen of each of serotypes 1 to 4, wherein each of said dengue antigens is a live attenuated chimeric dengue virus, preferably a chimeric YF/dengue virus, more preferably a chimeric YF/dengue virus which comprises an attenuated YF genomic backbone whose prM-E sequence has been substituted with the prM-E sequence of dengue virus.

[0057] Preferably, a live attenuated chimeric dengue virus to be analysed according to the present invention comprises one or more proteins from a dengue virus and one or more proteins from a different flavivirus. Preferably, the different

flavivirus is an yellow fever virus (i.e. a chimeric YF/dengue virus). Preferably a live attenuated chimeric dengue virus to be analysed according to the present invention comprises an attenuated yellow fever virus genome whose prM-E sequence has been substituted with the prM-E sequence of a dengue virus. Alternatively, a live attenuated chimeric dengue virus to be analysed according to the present invention comprises one or more proteins from a first dengue virus and one or more proteins from a second dengue virus (i.e. a chimeric dengue/dengue virus). Preferably said first dengue virus and said second dengue virus are of different serotypes. Where said first dengue virus and said second dengue virus are of the same serotype, said first and second dengue viruses are different strains.

[0058] If necessary, the virus concentration in the sample may be increased by growing the virus in suitable host cells. Suitable host cells include mammalian cells such as e.g. Vero cells. It is preferred that Vero cells are seeded in T25 cm² flasks at 1.5e⁶ cells/flask and grown overnight to 80-90% confluency by incubating at 37° C. with 5% CO₂ tension. Preferably, the virus is diluted in assay medium DMEM containing 1% FBS and 1% Penicillin/Streptomycin. Preferably, Vero cells are infected with TDV at 0.1 MOI and allowed virus adsorption for 1 hour at 37° C. 5% CO₂ incubator. For obtaining high quality viral RNA a preferred method is described in the Examples section.

[0059] "Purifying virus nucleic acids from the sample" may be performed using any method for DNA or RNA isolation known in the art. Any method known in the art may be used for DNA extraction or purification. Suitable methods comprise inter alia steps such as centrifugation steps, precipitation steps, chromatography steps, dialyzing steps, heating steps, cooling steps and/or denaturation steps.

[0060] "Preparing double-stranded DNA templates" herein means that if the virus nucleic acid is an RNA, the method further includes a step of reverse transcription with a reverse transcriptase, a reverse primer and a mixture of nucleotides such that a RNA/DNA hybrid is generated which is subsequently converted into a double-stranded cDNA template. The reaction may be carried out in an automated thermal cycler. Preferably, for RNA purified from Dengue virus the reverse transcription is performed using as reverse primer an oligonucleotide having the nucleotide sequence of SEQ ID NO:6 at a temperature in the range from about 45° C. to about 65° C., preferably at a temperature in the range from about 50° C. to about 55° C.

[0061] If the virus nucleic acid is a single-stranded DNA, the double-stranded DNA-template is generated by the activity of a template-dependent DNA polymerase. If the virus nucleic acid purified from the virus is already a double-stranded nucleic acid, there is no need for this optional step.

[0062] "Amplifying said viral nucleic acids of step a) or step b) such that PCR products having a length of at least 2 kb and comprising at least one attenuation locus are generated" herein means that the virus nucleic acids of a) or the double-stranded DNA templates of step b) are incubated with suitable primers (forward and reverse primers), nucleotides and a template-dependent DNA polymerase such that a PCR product of at least 2 kb is generated, wherein the PCR product comprises one or more of the attenuation loci of the virus. The reaction may be carried out in an automated thermal cycler.

[0063] The term “primer” as used herein denotes an oligonucleotide that acts as an initiation point of nucleotide synthesis under conditions in which synthesis of a primer extension product complementary to a nucleic acid strand is induced. The primers may be generated such that they are able to bind to the sequence to be reverse transcribed or amplified. The primer may comprise from 5 to 50 nucleotides, preferably, 8 to 30 nucleotides; more preferably the primer comprises 10 to 20 nucleotides. The primer sequences are complementary to at least a portion of the sequence to be reverse transcribed or amplified. Preferably, the primer sequences are fully complementary to a sequence within the sequence to be reverse transcribed or amplified. Methods for the design of sequence specific primers and probes are known in the art.

[0064] Preferably, the PCR product has a length of at least 3 kb, and more preferably 5 kb length. Preferably, the PCR product comprises two attenuation loci. More preferably, the PCR product comprises three attenuation loci. For TDV, the PCR product preferably comprises the positions 57, 2579 and/or 5270 of the TDV2 sequence (SEQ ID NO:23).

[0065] Preferably, the PCR amplification reaction of step c) comprises at least 10 PCR cycles using a DNA polymerase, suitable forward and reverse primers and nucleotides. More preferably, 10 to 50 cycles, even more preferred 10 to 30 cycles are used. For Dengue virus PCR product amplification, the PCR amplification reaction comprises the use of a first primer pair having the nucleotide sequence of SEQ ID NO: 9 (5F) and a second primer having the nucleotide sequence of SEQ ID NO: 6 (3R). It is also preferred that the PCR amplification reaction includes a denaturing step at a temperature from about 90° C. to about 99° C., an annealing step at a temperature from about 45° C. to about 75° C., and an extension step at a temperature from about 60° C. to about 75° C., more preferably denaturing is performed at about 95° C., annealing is performed at 68° C. with a 1° C. drop for each cycle over 10 cycles and extension is performed at 68° C.

[0066] Preferably, after the amplification step c) and prior to the sequencing step d) the method additionally comprises the steps of:

[0067] i) end repairing of the PCR products;

[0068] ii) optionally, purifying the end repaired PCR product; and

[0069] iii) connecting the two complementary strands of the PCR product prepared in step c) with a linking oligonucleotide, wherein the 3' end of one complementary strand is linked to the 5' end of the other complementary strand via the linking oligonucleotide in each double stranded nucleic acid fragment, whereby the linking oligonucleotide provides a single stranded portion of the resulting linked nucleic acid fragments.

[0070] End repairing of the PCR products is intended to generate blunt-end PCR products for the subsequent ligation of adaptors for SMRT sequencing. End repairing can be achieved by the use of suitable DNA polymerases such as the T4 DNA polymerase (Fast DNA End Repair Kit from Thermofischer) and nucleotides.

[0071] “Connecting the two complementary strands of the PCR product prepared in step c) with a linking oligonucleotide” herein means that by ligating of hairpin adaptors to both ends of the above double-stranded DNA product a closed, single-stranded circular DNA, the so-called SMRTbell, is created. Suitable hairpin adaptors are available from

Pacific Biosciences (SMRTbell® Template Prep Kit). The hairpin adaptor is linked to the blunt-ended PCR product using a suitable ligase.

[0072] Preferably, prior to step d) the amplified PCR products are purified. Methods for purifying of PCR products are known in the art. Suitable methods include the separation on an agarose gel or purification by column or beads.

[0073] “Single-molecule real-time (SMRT) sequencing” of step d) herein includes the sequencing techniques and platforms provided by the company Pacific Biosciences. SMRT sequencing preferably includes the PacBio RS and RSII platform and the PacBio Sequel platform, most preferably SMRT sequencing is performed using the PacBio RSII platform. When a sample with the SMRTbell structure is loaded to a chip called a SMRTcell, a SMRTbell diffuses into a sequencing unit called a zero-mode waveguide (ZMW), which provides the smallest available volume for light detection. The ZMW is a structure that creates an illuminated observation volume that is small enough to observe only a single nucleotide of DNA being incorporated by DNA polymerase.

[0074] For carrying out the SMRT sequencing to this SMRT cell a primer complementary to the hairpin adaptor is added with a mixture of mononucleotides under suitable conditions. Each of the four mononucleotide bases is attached to one of four different fluorescent dyes. When a nucleotide is incorporated by the DNA polymerase, the fluorescent tag is cleaved off and diffuses out of the observation area of the ZMW where its fluorescence is no longer observable. A detector detects the fluorescent signal of the nucleotide incorporation, and the base call, i.e. the correlation of the base with the fluorescent signal, is made according to the corresponding fluorescence of the dye.

[0075] SMRT sequencing herein also includes a technique developed by the company Pacific Biosciences. Single-molecule real-time (SMRT) sequencing is a parallelized single molecule DNA sequencing method. Single-molecule real-time sequencing utilizes a zero-mode waveguide (ZMW). A single DNA polymerase enzyme is affixed at the bottom of a ZMW with a single molecule of DNA as a template. The ZMW is a structure that creates an illuminated observation volume that is small enough to observe only a single nucleotide of DNA being incorporated by DNA polymerase. Each of the four DNA bases is attached to one of four different fluorescent dyes. When a nucleotide is incorporated by the DNA polymerase, the fluorescent tag is cleaved off and diffuses out of the observation area of the ZMW where its fluorescence is no longer observable. A detector detects the fluorescent signal of the nucleotide incorporation, and the base call is made according to the corresponding fluorescence of the dye; see Eid et al., *Science* 323 (2009), 133-138.

[0076] Preferably, SMRT sequencing is carried out using the PacBio RS and RSII device and the PacBio Sequel device, most preferably SMRT sequencing is performed using the PacBio RSII device. Pacific Biosciences has recently introduced a further improved version of the Sequel device, designated as Sequel II which can also be used for carrying out the present invention. For the data analysis, Pacific Biosciences offers SMRT analysis software which can be used for performance of the present invention.

[0077] Different chemistries are commercially available for the DNA bases attached to the fluorescent dyes. Orig-

nally, the so-called C1 chemistry was offered, followed by the so-called P6-C4 chemistry offered by Pacific Biosciences. All of them are encompassed by the scope of the invention. DNA polymerases have been optimized over time for use in the PacBio system. All of them are encompassed by the invention.

[0078] In order to reduce the sequence error rate of SMRT sequencing, a preferably high number of sequence reads is used for data analysis. More preferably, at least 10^3 , more preferably, 10^4 , most preferred 5×10^4 reads/run are used. The higher the number of reads confirming a particular sequence, the higher is the statistical reliability of this given sequence. Reads below 10 are considered as unreliable.

[0079] SMRT sequencing differs from second-generation sequencing (SGS) technologies such as those provided by the companies Illumina (Miniseq, MiSeq, NextSeq or HiSeq) or Ion Torrent (PGM, S5, Proton), and Roche (454, pyrosequencing). The SGS technologies are characterized by much shorter read length of up to 400 bp, while having a much lower sequence error rate than SMRT sequencing. Quantifying low-frequency revertants in oral poliovirus vaccine was performed using SGS technology, in particular MiSeq System (Illumina) or HiSeq 1500 System (Illumina), and the accuracy was shown to be comparable with the prior art MAPREC (mutant analysis by polymerase chain reaction and restriction enzyme cleavage) assay (Sarcey et al., J. Virol. Methods 246 (2017), 75-80). In view of the short read-lengths SGS technologies do not provide an individual sequence of a large insert such as 1 kb or more, but rather for such inserts the SGS technologies provide an average sequence of different polynucleotides.

[0080] The wild-type virus sequences for different viruses are available from public databases are e.g. hosted by NCBI such as Entrez (www.ncbi.nlm.nih.gov). The wild-type virus sequence of Dengue virus serotype 2 differs from the TDV2 nucleotide sequence set forth in SEQ ID NO:23 in that the wild-type sequence has nucleotide C at position 57, nucleotide G at position 2579 and nucleotide A on position 5270.

[0081] “Comparing at least one nucleotide at said at least one attenuation locus in said PCR products with the corresponding nucleotide in the wild-type virus sequence” herein means a comparison of the obtained sequence of the PCR products with the wild-type sequence at the attenuation locus/loci. The comparison can be performed manually or using a computer program. If the obtained sequence of the PCR products has the same nucleotide at the attenuation locus as the wild-type, this is indicative of the presence of a revertant in the sample. If the obtained sequence of the PCR products has the same nucleotide at the attenuation locus as the attenuated virus, this is indicative that no revertant virus is present in said sample.

[0082] In one embodiment, if the attenuated virus is TDV, at least one nucleotide at said at least one attenuation locus in said PCR products is compared with the corresponding nucleotide in the wild-type virus sequence, preferably at least one nucleotide at two attenuation loci in said PCR products is compared with the corresponding nucleotide in the wild-type virus sequence and more preferably at least one nucleotide at three attenuation loci in said PCR products is compared with the corresponding nucleotide in the wild-type virus sequence.

[0083] In one embodiment, if the attenuated virus is TDV, at least one nucleotide at at least one attenuation locus selected from nucleotides corresponding to nucleotides 57,

2579 and 5270 in the nucleic acid sequence according to SEQ ID No. 23 in said PCR products is compared with the corresponding nucleotide in the wild-type virus sequence. Preferably at least one nucleotide at two attenuation loci selected from nucleotides corresponding to nucleotides 57, 2579 and 5270 in the nucleic acid sequence according to SEQ ID No. 23 in said PCR products is compared with the corresponding nucleotide in the wild-type virus sequence.

[0084] More preferably at least one nucleotide at three attenuation loci at from nucleotides corresponding to nucleotides 57, 2579 and 5270 in the nucleic acid sequence according to SEQ ID No. 23 in said PCR products is compared with the corresponding nucleotide in the wild-type virus sequence.

[0085] According to a further aspect of the invention is provided a method for the quantitative analysis of a mixture of virus haplotypes in a virus sample comprising the steps:

[0086] a) purifying virus nucleic acids from said sample;

[0087] b) preparing double-stranded DNA templates from said virus nucleic acids by using a primer comprising a unique molecular identifier (UMI) sequence such that said double-stranded DNA templates comprise the UMI sequence;

[0088] c) amplifying said double-stranded DNA templates by polymerase chain reaction (PCR) such that PCR products having a length of at least 2 kb and comprising a unique molecular identifier (UMI) sequence are generated;

[0089] d) sequencing said PCR products by single-molecule real-time (SMRT) sequencing; and

[0090] e) determining the relative amounts of said PCR products by analyzing the UMI count and the Read count thereby the virus haplotypes contained in the sample.

[0091] The embodiments for the virus and the sample are as defined above. A “virus haplotype” relates to the entire nucleotide sequence of a given virus. A virus haplotype may be defined as the combination of nucleotides at particular loci such as the attenuation loci. For example, the Dengue virus haplotype T:57-A:2579-T:5270 is different from the haplotype T:57-A:2579-C:5270 by the difference in nucleotide 5270.

[0092] “Unique molecular identifiers (UMIs)” are sequences of nucleotides applied to or identified in DNA molecules that may be used to distinguish individual DNA molecules from one another (see, e.g., Kivioja, Nature Methods 9, 72-74 (2012)). UMIs may be sequenced along with the DNA molecules with which they are associated to determine whether the read sequences are those of one source DNA molecule or another. The term “UMI” is used herein to refer to both the sequence information of a polynucleotide and the physical polynucleotide per se. The term “UMI” includes bar codes which are commonly used to distinguish reads of one sample from reads of other samples.

[0093] Commonly, multiple copies of a single source molecule are sequenced.

[0094] The source molecule for the sequencing may be PCR amplified before delivery to a flow cell. Each molecule in a PCR-amplified cluster is derived from the same source DNA molecule, but is separately sequenced. For error correction and other purposes, it can be important to determine that all reads from a single cluster are identified as deriving from the same source molecule. UMIs allow this grouping.

A DNA molecule that is copied by amplification or otherwise to produce multiple copies of the DNA molecule is referred to as a source DNA molecule.

[0095] UMIs may be applied to or identified in individual DNA molecules. In some embodiments, the UMIs may be applied to the DNA molecules by methods that physically link or bind the UMIs to the DNA molecules, e.g., by ligation or transposition through polymerase, endonuclease, transposases, etc. These “applied” UMIs are therefore also referred to as physical UMIs. In some contexts, they may also be referred to as exogenous UMIs. The UMIs identified within source DNA molecules are referred to as virtual UMIs. In some context, virtual UMIs may also be referred to as endogenous UMIs.

[0096] Physical UMIs may be defined in many ways. For example, they may be random, pseudo-random or partially random, or non-random nucleotide sequences that are inserted in adapters or otherwise incorporated in source DNA molecules to be sequenced. In some embodiments, the physical UMIs may be so unique that each of them is expected to uniquely identify any given source DNA molecule present in a sample. The collection of adapters is generated, each having a physical UMI, and those adapters are attached to fragments or other source DNA molecules to be sequenced, and each individual sequenced molecule has a UMI that helps to distinguish it from all other fragments. In such embodiments, a very large number of different physical UMIs (e.g., many thousands to millions) may be used to uniquely identify DNA fragments in a sample.

[0097] The physical UMI must have a sufficient length to ensure this uniqueness for each and every source DNA molecule. In some embodiments, a less unique molecular identifier can be used in conjunction with other identification techniques to ensure that each source DNA molecule is uniquely identified during the sequencing process. In such embodiments, multiple fragments or adapters may have the same physical UMI. Other information such as alignment location or virtual UMIs may be combined with the physical UMI to uniquely identify reads as being derived from a single source DNA molecule/fragment. In some embodiments, adapters include physical UMIs limited to a relatively small number of non-random sequences, e.g., 96 non-random sequences. Such physical UMIs are also referred to as non-random UMIs. In some embodiments, the non-random UMIs may be combined with sequence position information and/or virtual UMIs to identify reads attributable to a same source DNA molecule. The identified reads may be collapsed to obtain a consensus sequence that reflects the sequence of the source DNA molecule as described herein.

[0098] A “random UMI” may be considered a physical UMI selected as a random sample, from a set of UMIs consisting of all possible different oligonucleotide sequences given one or more sequence lengths. For instance, if each UMI in the set of UMIs has n nucleotides, then the set includes 4^n UMIs having sequences that are different from each other. A random sample selected from the 4^n UMIs constitutes a random UMI.

[0099] Conversely, a “non-random UMI” as used herein refers to a physical UMI that is not a random UMI. In some embodiments, available non-random UMIs are predefined for a particular experiment or application. In certain embodiments, rules are used to generate sequences to obtain a non-random UMI. For instance, the sequences of a set may

be generated such that the sequences have a particular pattern or patterns. In some embodiments, each sequence differs from every other sequence in the set by a particular number of (e.g., 2, 3, or 4) nucleotides. That is, no non-random UMI sequence can be converted to any other available non-random UMI sequence by replacing fewer than the particular number of nucleotides. In some embodiments, a non-random UMI is selected from a set of UMIs including fewer than all possible UMIs given a particular sequence length. For instance, a nonrandom UMI having 6 nucleotides may be selected from a total of 96 different sequences (instead of a total of $4^6=4096$ possible different sequences). In other embodiments, sequences are not randomly selected from a set. Instead, some sequences are selected with higher probability than other sequences.

[0100] The attachment of UMIs can be performed by adding a UMI sequence such as a random sequence to the 5' end of a forward primer or the 3' end of a reverse primer to be used in the PCR amplification reaction. Thereby, the UMI sequence is incorporated into the PCR amplification product.

[0101] If the virus nucleic acids are isolated virus RNA, step b) comprises a step of reverse transcription of the isolated virus RNA in the presence of a reverse transcriptase using a reverse primer additionally comprising at its 5' end a UMI sequence. Preferably, in case of Dengue virus, reverse transcription is performed using as reverse primer an oligonucleotide having the nucleotide sequence of SEQ ID NO:22 (UMI 3R) at a temperature in the range from about 45° C. to about 65° C., preferably at a temperature in the range from about 50° C. to about 55° C.

[0102] Alternatively, the UMI sequence can be incorporated into the forward or reverse primer at the PCR amplification of step c).

[0103] Preferably, the PCR amplification of step c) includes a PCR amplification reaction comprising at least 10 PCR cycles using a DNA polymerase and forward and reverse primers. More preferably, 10 to 50 PCR cycles, even more preferred 10 to 30 PCR cycles are performed. In case of Dengue virus, the PCR amplification reaction comprises the use of a first primer having the nucleotide sequence of SEQ ID NO:9 (5F) and a second primer having the nucleotide sequence of SEQ ID NO:21 (Primer IIA).

[0104] Preferably, the PCR amplification reaction of step c) includes a denaturing step at a temperature from about 90° C. to about 99° C., an annealing step at a temperature from about 45° C. to about 75° C., and an extension step at a temperature from about 60° C. to about 75° C.

[0105] Preferably, after the amplification step c) and prior to the sequencing step d) the method additionally comprises the steps of:

[0106] i) end repairing of the PCR products;

[0107] ii) optionally, purifying the end repaired PCR product; and

[0108] iii) connecting the two complementary strands of the PCR product prepared in step c) with a linking oligonucleotide, wherein the 3' end of one complementary strand is linked to the 5' end of the other complementary strand via the linking oligonucleotide in each double stranded nucleic acid fragment, whereby the linking oligonucleotide provides a single stranded portion of the resulting linked nucleic acid fragments.

[0109] Embodiments of steps c) and d) are as defined above.

[0110] “Read count” represents the number of reads. It is a category of quantification scheme commonly employed in single molecule sequencing. Although the read count-based scheme is similar to the common approaches used for bulk RNA-sequencing, the miniscule quantity of transcripts captured from a single cell requires cDNA amplification for library construction; this inevitably results in large amplification bias. Thus, if a truncated or erroneous cDNA fragment is generated which binds stronger to the used primer pair, the truncated or erroneous cDNA fragment will be amplified in a higher amount compared to the correct DNA fragment. “UMI count” relates to number of the uniquely mapped reads only. To mitigate the bias in read counts, several recent scRNA-sequencing protocols have employed an additional step in which individual transcripts are modified with unique molecular identifiers (UMIs) before amplification, resulting in a more accurate quantification of the transcript count. The UMI count is determined by considering only the Read count of sequences which contain a unique UMI sequence.

[0111] According to a preferred embodiment the method of the present invention comprises prior to step d) the steps:

[0112] i) determining the sequence of the PCR products at at least one attenuation locus separately by Sanger sequencing;

[0113] ii) comparing at least one nucleotide at said at least one attenuation locus with the wild-type virus sequence; and

[0114] iii) if the wild-type sequence is determined at at least one attenuation locus, sequencing said PCR products by SMRT sequencing of step d) is carried out.

[0115] Sanger sequencing, i.e. a method of DNA sequencing based on the selective incorporation of chain-terminating dideoxynucleotides by DNA polymerase during in vitro DNA replication is significantly less expensive compared to SMRT sequencing. If the Sanger sequencing at the individual attenuation loci indicates the presence of a revertant virus, SMRT sequencing of the entire virus nucleic acid molecule is conducted to discriminate single, double and triple revertants at the attenuation loci.

[0116] A suitable primer pair for Sanger sequencing of the first attenuation locus at position 57 of the TDV2 nucleotide sequence set forth as SEQ ID NO:23 is the forward primer having the nucleotide sequence of SEQ ID NO: 13 and the reverse primer of SEQ ID NO: 17. A suitable primer pair for Sanger sequencing of the second attenuation locus at position 2579 is the forward primer having the nucleotide sequence of SEQ ID NO: 14 and the reverse primer of SEQ ID NO: 18. Suitable primer pairs for Sanger sequencing of the third attenuation locus at position 5270 are forward primers having a nucleotide sequence of SEQ ID Nos: 15 or 16, and as reverse primers having a nucleotide sequence of SEQ ID Nos: 19 or 20.

[0117] According to a further aspect the present invention provides the use of the method according to the first aspect in the quality control of vaccines containing live, attenuated virus. For quality control of vaccine samples, it is of particular importance to detect the existence of revertants in the vaccine sample. It is readily apparent that dengue virus having a reversion at each of the three attenuation loci (triple revertant) has the wild-type replication capacity and may therefore be virulent. However, a vaccine composition com-

prising a double revertant having a reversion at two of the attenuation loci and a monorevertant having a reversion at the third attenuation locus may also pose a risk in view of the possibility of recombination events between the double revertant and the single revertant which may lead to increased virulence.

[0118] According to a further aspect the present invention provides the use of the method according to the second aspect in the diagnosis of a virus infection in a patient sample. The method according to the second aspect may be used for the diagnosis of a virus infection in a virus infected patient. The method also enables the detection of co-infection with different subtypes of viruses such as different dengue serotypes. This may provide advantages in the personalized treatment of the infected patients. Preferably, the method for the diagnosis of a virus infection in a virus infected patient is carried out in vitro on a biological sample removed from the patient. Preferably, the sample is a blood sample, more preferably the sample is a serum sample.

[0119] According to a further aspect of the present invention a kit is provided for detecting the presence of at least one revertant virus in a sample containing an attenuated dengue virus comprising a first primer selected from an oligonucleotide having the nucleotide sequence of SEQ ID NO: 1, SEQ ID NO: 3, SEQ ID NO: 5, SEQ ID NO: 7, SEQ ID NO:9 and SEQ ID NO:11, and comprising a second primer selected from an oligonucleotide having the nucleotide sequence of SEQ ID NO: 2, SEQ ID NO: 4, SEQ ID NO: 6, SEQ ID NO:8, SEQ ID NO: 10 and SEQ ID NO: 12. Preferably, the first primer comprises the nucleotide sequence of SEQ ID NO:9 (5F) and the second primer comprises the nucleotide sequence of SEQ ID NO:6 (3R). More preferably, the first primer consists of the nucleotide sequence of SEQ ID NO:9 and the second primer consists of the nucleotide sequence of SEQ ID NO:6.

[0120] The kit may further include a suitable reverse transcriptase, nucleotides and suitable buffers known in the art.

[0121] According to a further aspect of the present invention a kit is provided for the quantitative analysis of a mixture of virus haplotypes in a dengue virus sample comprising a first primer comprising the nucleotide sequence of SEQ ID NO:9 (5F) and a second primer comprising the nucleotide sequence of SEQ ID NO:21 (Primer IIA). Preferably the first primer consists of the nucleotide sequence of SEQ ID NO:9 and the second primer consists of the nucleotide sequence of SEQ ID NO:21.

[0122] The kit may further include a suitable reverse transcriptase, nucleotides and suitable buffers known in the art.

[0123] The invention is further described in the following examples which are solely for the purpose of illustrating specific embodiments of the invention, and are also not to be construed as limiting the scope of the invention in any way.

LIST OF EMBODIMENTS

[0124] 1. A method for detecting the presence of at least one revertant virus in a sample containing an attenuated virus comprising the steps:

[0125] a) purifying virus nucleic acids from said sample;

[0126] b) optionally, preparing double-stranded DNA templates from said virus nucleic acids;

- [0127] c) amplifying said viral nucleic acids of a) or said double-stranded DNA templates of b) by polymerase chain reaction (PCR) such that PCR products having a length of at least 2 kb and comprising at least one attenuation locus of said attenuated virus are generated;
- [0128] d) sequencing said PCR products by single-molecule real-time (SMRT) sequencing; and
- [0129] e) comparing at least one nucleotide at said at least one attenuation locus in said PCR products with the corresponding nucleotide in the wild-type virus sequence, wherein the presence of the wild-type nucleotide at said at least one attenuation locus is indicative of the presence of a revertant virus in said sample.
- [0130] 2. The method of embodiment 1, wherein said virus is an RNA virus, preferably wherein the RNA virus is selected from dengue virus, poliovirus, rubella virus, measles virus, yellow fever virus, mumps virus and influenza virus, more preferably the RNA virus is dengue virus.
- [0131] 3. The method of embodiment 1 or 2, wherein the at least one revertant virus is present in the sample containing the attenuated virus in a quantity of less than 1%, preferably less than 0.2% based on the total amount of virus in the sample.
- [0132] 4. The method of any one of embodiments 1 to 3, wherein the at least one revertant virus differs from the attenuated virus by less than 10 nucleotides, preferably less than 5 nucleotides, more preferably by less than 3 nucleotides, and most preferably by one nucleotide.
- [0133] 5. The method of embodiment 1, wherein said virus is a DNA virus, preferably wherein the DNA virus is varicella zoster virus, vaccinia virus, smallpox virus, Herpes simplex virus and chikungunya virus.
- [0134] 6. The method of any one of embodiments 1 to 5, wherein said PCR products have a length of at least 3 kb, preferably of at least 5 kb.
- [0135] 7. The method of any one of embodiments 1 to 6, wherein said PCR products comprise at least two attenuation loci, preferably at least three attenuation loci.
- [0136] 8. The method of any one of embodiments 1, 2, 6 or 7, wherein said attenuated virus is a live, attenuated dengue virus and said at least one attenuation locus is selected from positions 57, 2579 and/or 5270 of the nucleotide sequence of dengue virus type 2 according to SEQ ID NO: 23.
- [0137] 9. The method of any one of embodiments 1 to 8, wherein prior to step a) the live, attenuated virus is grown in Vero cells.
- [0138] 10. The method of any one of embodiments 1, 2 or 6 to 9, wherein the virus nucleic acids are isolated virus RNA and step b) comprises a step of reverse transcription of the isolated virus RNA in the presence of a reverse transcriptase using a reverse primer.
- [0139] 11. The method of embodiment 10, wherein the reverse transcription is performed using as reverse primer an oligonucleotide comprising the nucleotide sequence of SEQ ID NO:6 (3R) at a temperature in the range from about 45° C. to about 65° C., preferably at a temperature in the range from about 50° C. to about 55° C., preferably the oligonucleotide consists of the nucleotide sequence of SEQ ID NO:6.
- [0140] 12. The method of any one of embodiments 1 to 11, wherein step c) includes a PCR amplification reaction comprising at least 10 PCR cycles using a DNA polymerase and forward and reverse primers.
- [0141] 13. The method of embodiment 12, wherein the PCR amplification reaction comprises the use of a first primer comprising the nucleotide sequence of SEQ ID NO:6 (3R) and a second primer comprising the nucleotide sequence of SEQ ID NO:9 (5F), preferably the first primer consists of the nucleotide sequence of SEQ ID NO:6 and the second primer consists of the nucleotide sequence of SEQ ID NO:9.
- [0142] 14. The method of embodiment 12 or 13, wherein the PCR amplification reaction includes a denaturing step at a temperature from about 90° C. to about 99° C., an annealing step at a temperature from about 45° C. to about 75° C., and an extension step at a temperature from about 60° C. to about 75° C.
- [0143] 15. The method of any one of embodiments 1 to 14, wherein the amplified PCR products obtained in step c) are purified prior to step d).
- [0144] 16. The method of any one of embodiments 1 to 15, wherein prior to step d) the amplified PCR products are provided as a library suitable for SMRT sequencing.
- [0145] 17. The method of embodiment 16, wherein prior to step d) the method additionally comprises:
- [0146] i) end repairing of the PCR products;
- [0147] ii) optionally, purifying the end repaired PCR product; and
- [0148] iii) connecting the two complementary strands of the PCR product prepared in step c) with a linking oligonucleotide, wherein the 3' end of one complementary strand is linked to the 5' end of the other complementary strand via the linking oligonucleotide in each double stranded nucleic acid fragment, whereby the linking oligonucleotide provides a single stranded portion of the resulting linked nucleic acid fragments.
- [0149] 18. A method for the quantitative analysis of a mixture of virus haplotypes in a virus sample comprising the steps of
- [0150] a) purifying virus nucleic acids from said sample;
- [0151] b) preparing double-stranded DNA templates from said virus nucleic acids by using a primer comprising a unique molecular identifier (UMI) sequence such that said double-stranded DNA templates comprise the UMI sequence;
- [0152] c) amplifying said double-stranded DNA templates by polymerase chain reaction (PCR) such that PCR products having a length of at least 2 kb and comprising a unique molecular identifier (UMI) sequence are generated;
- [0153] d) sequencing said PCR products by single-molecule real-time (SMRT) sequencing; and
- [0154] e) determining the relative amounts of said PCR products by analysing the UMI count and the Read count, thereby quantifying the virus haplotypes contained in the sample.
- [0155] 19. The method of embodiment 18, wherein any one of the virus haplotypes is present in the virus sample in a quantity of less than 10%, preferably less

- than 1%, more preferably less than 0.2% based on the total amount of virus in the sample.
- [0156] 20. The method of embodiment 18, wherein the virus is an RNA virus, preferably wherein the RNA virus is selected from dengue virus, poliovirus, rubella virus, measles virus, yellow fever virus, mumps virus and influenza virus, more preferably the RNA virus is dengue virus.
- [0157] 21. The method of any one of embodiments 18 to 20, wherein at least two of the virus haplotypes differ from each other by less than 10 nucleotides, preferably less than 5 nucleotides, more preferably by less than 3 nucleotides, and most preferably by one nucleotide.
- [0158] 22. The method of embodiment 18, wherein the virus is a DNA virus, preferably wherein the DNA virus is varicella zoster virus.
- [0159] 23. The method of any one of embodiments 18 to 21, wherein said virus is a live, attenuated dengue virus, preferably a dengue virus comprising at least one attenuation locus selected from positions 57, 2579 and/or 5270 of the nucleotide sequence of dengue virus type 2 according to SEQ ID NO: 23.
- [0160] 24. The method of any one of embodiments 18 to 23, wherein said PCR products have a length of at least 3 kb, preferably of at least 5 kb.
- [0161] 25. The method of any one of embodiments 18, 19 or 23, wherein the virus nucleic acids are isolated virus RNA and step b) comprises a step of reverse transcription of the isolated virus RNA in the presence of a reverse transcriptase using a reverse primer additionally comprising at its 5' end a UMI sequence.
- [0162] 26. The method of embodiment 25, wherein the reverse transcription is performed using as reverse primer an oligonucleotide comprising the nucleotide sequence of SEQ ID NO:22 (UMI 3R) at a temperature in the range from about 45° C. to about 65° C., preferably the oligonucleotide consists of the nucleotide sequence of SEQ ID NO:22 and the temperature is in the range from about 50° C. to about 55° C.
- [0163] 27. The method of any one of embodiments 18 to 26, wherein step c) includes a PCR amplification reaction comprising at least 10 PCR cycles using a DNA polymerase and forward and reverse primers.
- [0164] 28. The method of embodiment 27, wherein the PCR amplification reaction comprises the use of a first primer comprising the nucleotide sequence of SEQ ID NO:21 (Primer IIA) and a second primer comprising the nucleotide sequence of SEQ ID NO: 9 (5F), preferably the first primer consists of the nucleotide sequence of SEQ ID NO:21 and the second primer consists of the nucleotide sequence of SEQ ID NO:9.
- [0165] 29. The method of embodiment 27 or 28, wherein the PCR amplification reaction includes a denaturing step at a temperature from about 90° C. to about 99° C., an annealing step at a temperature from about 45° C. to about 75° C., and an extension step at a temperature from about 60° C. to about 75° C.
- [0166] 30. The method of any one of embodiments 18 to 29, wherein the amplified PCR products obtained in step c) are purified prior to step d).
- [0167] 31. The method of any one of embodiments 18 to 30, wherein prior to step d) the amplified PCR products are provided as a library suitable for SMRT sequencing.
- [0168] 32. The method of embodiment 31, wherein after the amplification step c) and prior to the sequencing step d) the method additionally comprises:
- [0169] i) end repairing of the PCR products;
- [0170] ii) optionally, purifying the end repaired PCR product; and
- [0171] iii) connecting the two complementary strands of the double stranded nucleic acid fragments prepared in step (b) with a linking oligonucleotide, wherein the 3' end of one complementary strand is linked to the 5' end of the other complementary strand via the linking oligonucleotide in each double stranded nucleic acid fragment, whereby the linking oligonucleotide provides a single stranded portion of the resulting linked nucleic acid fragments.
- [0172] 33. The method according to any one of embodiments 1 to 32, wherein said sample is a composition comprising an attenuated virus, or a sample from a patient.
- [0173] 34. The method of any one of embodiments 1 to 33, wherein prior to step d) the method comprises:
- [0174] i) determining the sequence of the PCR products at at least one of the attenuation locus separately by Sanger sequencing;
- [0175] ii) comparing at least one nucleotide at said at least one attenuation locus with the wild-type virus sequence; and
- [0176] iii) if the wild-type sequence is determined at at least one attenuation locus, sequencing said PCR products by SMRT sequencing of step d) is carried out.
- [0177] 35. Use of the method according to any one of embodiments 1 to 17 in the quality control of vaccines containing live, attenuated virus.
- [0178] 36. Use of the method according to any one of embodiments 18 to 34 in the diagnosis of a virus infection in a patient sample.
- [0179] 37. A kit for detecting the presence of at least one revertant virus in a sample containing an attenuated dengue virus comprising a first primer selected from an oligonucleotide having the nucleotide sequence of SEQ ID NO: 1, SEQ ID NO: 3, SEQ ID NO: 5, SEQ ID NO: 7, SEQ ID NO:9 and SEQ ID NO:11, and comprising a second primer selected from an oligonucleotide having the nucleotide sequence of SEQ ID NO: 2, SEQ ID NO: 4, SEQ ID NO: 6, SEQ ID NO:8, SEQ ID NO: 10 and SEQ ID NO: 12, preferably the first primer comprises the nucleotide sequence of SEQ ID NO: 9 (5F) and the second primer comprises the nucleotide sequence of SEQ ID NO: 6 (3R), more preferably, the first primer consists of the nucleotide sequence of SEQ ID NO:9 and the second primer consists of the nucleotide sequence of SEQ ID NO:6.
- [0180] 38. A kit for the quantitative analysis of a mixture of virus haplotypes in a dengue virus sample comprising a first primer comprising the nucleotide sequence of SEQ ID NO:21 (Primer IIA) and a second primer comprising the nucleotide sequence of SEQ ID NO: 9 (5F), preferably the first primer consists of the nucleotide sequence of SEQ ID NO:21 and the second primer consists of the nucleotide sequence of SEQ ID NO:9.

EXAMPLES

Material and Methods

Growth and Infection of Vero Cells

[0181] Vero WHO cells were seeded in T25 cm² flasks at 1.5e⁶ cells/flask and grown overnight to 80-90% confluency by incubating at 37° C. with 5% CO₂ tension. Cells were observed for health and confluency under an inverted microscope before the start of the infection experiment.

[0182] TDV-2 virus was diluted in assay medium DMEM containing 1% FBS and 1% Penicillin/Streptomycin (HyClone, Cat #SV30010). Vero cells were infected with TDV-2 at 0.1 MOI and allowed virus adsorption for 1 hour at 37° C. 5% CO₂ incubator. Upon completion of incubation, inoculum was aspirated and cells were replenished with 6 ml of fresh assay medium. The flasks were allowed to continue incubate at this temperature for additional 5-7 days. Medium was changed on day 5 and virus was harvested on day 7 by collecting and spinning the cell culture supernatant at 10,000×g for 10 minutes in a 4° C. cooled centrifuge.

performed with selected primers followed by PCR. The PCR conditions were then optimized. The PCR products were enriched by nested PCR and further amplifications. The PCR products were treated with a kit for the production of pacbio library. Smaller products were removed using Ampure beads. Then SMRT sequencing was performed and data was analyzed using bioinformatic tools.

Primer Selection

[0185] The regions surrounding the variants at positions 57 and 5270 were analyzed using the RNA folding software (<http://unafold.ma.albany.edu/?q=mfold>) to identify regions having minimal secondary structure. Segments were chosen having unpaired regions and six primer pairs were chosen in these segments using the PRIMER3 server at NCBI. Since only 57 bp are available at the 5' end (FIG. 2) some of these primers overlap with one another (Table 1).

List of Primers Used for RT-PCR Optimization

TABLE 1

Primer sequences and characteristics. Note that the reverse primer is also the primer for reverse transcription. The primers are in pairs.						
Name	Dir.	Sequence	Start	Stop	Tm	Product
1F	F	AAGACAGATTCTTTGAGGGAGCTA	26	49	59	5368
1R	R	GGTCTGTGAAATGGGCTTCG	5393	5374	59	5368
2F	F	GACAGATTCTTTGAGGGAGCTAA	28	50	58	5367
2R	R	GGGTCTGTGAAATGGGCTTC	5394	5375	59	5367
3F	F	ACAGATTCTTTGAGGGAGCTAAG	29	51	58	5364
3R	R	GTCTGTGAAATGGGCTTCGTC	5392	5372	60	5364
4F	F	AGACAGATTCTTTGAGGGAGCTAAG	27	51	60	5366
4R	R	GTCTGTGAAATGGGCTTCGT	5392	5373	58	5366
5F	F	GACAGATTCTTTGAGGGAGCTAAG	28	51	59	5373
5R	R	CTTGCTGGGTCTGTGAAATGG	5400	5380	59	5373
6F	F	AAGACAGATTCTTTGAGGGAGCTAA	26	50	60	5370
6R	R	TGGGTCTGTGAAATGGGCTTC	5395	5375	61	5370

Viral RNA Extraction

[0183] 140 µl of TDV2/Vero WHO viral cell culture supernatant was subjected to RNA extraction using 3 different commercial kits. The extraction was carried out according to the instructions came with the kits. Qiagen QIAamp Viral RNA Mini Kit (Lot #145051122), Omega Viral RNA extraction kit and Roche's High Pure RNA Isolation kit as per the instructions of the manufacturers without addition of carrier RNA. The purified RNA yield was measured on a Nanodrop spectrophotometer.

Optimization of Pacbio Sequencing for TDV Analysis

Process Summary: Flow Chart

[0184] The pacbio sequencing work flow has been described in FIG. 1. Reverse transcription (RT) was first

TABLE 2

Correlation of Primer Sequences with SEQ ID numbers. Primers were obtained from IDT DNA (www.idtdna.com). The reverse primer in each case was used as the primer for reverse transcriptase in the first strand cDNA reaction. The pair was then used for PCR.	
Sequence	Name
AAGACAGATTCTTTGAGGGAGCTA	1F (SEQ ID NO: 1)
GGTCTGTGAAATGGGCTTCG	1R (SEQ ID NO: 2)
GACAGATTCTTTGAGGGAGCTAA	2F (SEQ ID NO: 3)
GGGTCTGTGAAATGGGCTTC	2R (SEQ ID NO: 4)

TABLE 2-continued

Correlation of Primer Sequences with SEQ ID numbers. Primers were obtained from IDT DNA (www.idtdna.com). The reverse primer in each case was used as the primer for reverse transcriptase in the first strand cDNA reaction. The pair was then used for PCR.	
Sequence	Name
ACAGATTCTTTGAGGGAGCTAAG	3F (SEQ ID NO: 5)
GTCTGTGAAATGGGCTTCGTC	3R (SEQ ID NO: 6)
AGACAGATTCTTTGAGGGAGCTAAG	4F (SEQ ID NO: 7)
GTCTGTGAAATGGGCTTCGT	4R (SEQ ID NO: 8)
GACAGATTCTTTGAGGGAGCTAAG	5F (SEQ ID NO: 9)
CTTGCTGGGTCTGTGAAATGG	5R (SEQ ID NO: 10)
AAGACAGATTCTTTGAGGGAGCTAA	6F (SEQ ID NO: 11)
TGGGTCTGTGAAATGGGCTTC	6R (SEQ ID NO: 12)

RT-PCR Optimization

[0186] Reverse transcription was carried out on 6.5 μL (139.1 ng) of the TDV2 viral RNA with the reverse primers (Table 1). The six primer pairs were used in touchdown PCR reactions. Touchdown PCR was done as follows. 95° C. for 30 sec., annealing at 66° C. for 20 sec. with a 1 degree drop each cycle for 10 cycles (i.e., cycle 1—66° C., cycle 2—65° C., cycle 3—64° C. etc.) and extension at 68° C. for 6 min. An additional 25 cycles were run at 95° C. for 30 sec., annealing at 53° C. for 20 sec. and extension at 68° C. for 6 minutes. RT-PCR products were run on 1% agarose gels in 1×TAE buffer.

[0187] The two best-performing primer pairs (1 and 4) were then used with nested PCR primers (5F and 3R, from the same set in Table 1) and optimized for number of cycles. Primer pair 1 (1F/R) with fewer non-specific bands was taken forward to annealing temperature optimization. Five temperatures were tested over a range of PCR cycles. The optimum of 11 cycles was chosen and a 0.45× Ampure XP bead purification was used to remove non-specific amplification products less than 3 kilobases.

PacBio Sequencing and Variant Analysis

[0188] The material from optimized RT-PCR reaction was used to construct a PacBio sequencing library. One SMRT-cell was sequenced, producing 366 Mbp of data and 12,467 reads of insert. The reads were aligned to the dengue virus reference genome (GenBank accession U87411) using BLASR and variants were called with QUIVER. The data are summarized in the Result section.

Identification of TDV2 Revertants

Sanger Sequencing to Confirm the Identity of the Revertants Four revertants (see Table 3 below) were amplified in Vero WHO cells and viral RNA was extracted using Omega Kit. Sanger sequencing was conducted to confirm identity of these TDV revertants.

TABLE 3

List of TDV2 revertants				
Revertant	Date produced	Reversion Site	TDV	WT
✓ D2/IC P1	Aug. 12, 1998	NS1-2579	A	G
✓ D2/IC P3	Jun. 9, 2004	NS3-5270	T	A
✓ D2/IC P5	Aug. 12, 1998	5'NC-57	T	C
✓ D2/IC P51	Aug. 12, 1998	5'NC-57	T	C
		NS1-2579	A	G

RT-PCR Optimization

[0189] Titan One tube RT-PCR kit (Roche) along with multiple sets of primers (not all shown here) and TDV2 viral RNA was used to optimize RT and PCR conditions to produce a single large product of around 5.3 kb size covering three attenuation loci including nt57, nt2579 and nt5270. Finally, a best pair of primers (shown below) that produced no non-specific products was selected.

D2-1 (Forward): (SEQ ID NO: 13)
AGTTGTTAGTCTACGTGGACCGAC

cD2-5358 (Reverse): (SEQ ID NO: 20)
GAAATGGGCTTCGTCCATGATAATCAGG

[0190] Optimized RT-PCR conditions were applied to four TDV2 revertants to amplify 5.3 kb product using 100 ng viral RNA per reaction tube. The RT-PCR was done as follows: 50° C. for 30 sec., 94° C. for 2 minutes (1 cycle), denaturation at 94° C. for 15 sec., annealing at 60° C. for 30 sec. and extension at 68° C. for 5 min (30 cycles) and final extension at 68° C. for 5 minutes (1 cycle). We successfully produced 5.3 kb long products that were subjected to sanger sequencing to identify TDV2 revertants.

Purification of PCR Product

[0191] RT-PCR reaction was optimized and PCR products of 5.35 kb were successfully amplified. Two purification kits from Qiagen were used to purify PCR products; Gel extraction and PCR purification kits. Gel extracted kit provided cleaner single band product and was selected for purification of PCR products.

Sanger Sequencing and Analysis

[0192] Purified PCR products were examined by Sanger sequencing and sequencing data files were edited and aligned to TDV2 reference sequencing at Takeda. The primers used in Sanger sequencing reactions and the positions of these primers on TDV2 genome covering three attenuation loci of interest can be seen in table and figure of Table 4 and FIG. 3, respectively.

List of Primers Used for Sanger Sequencing

Name	Direction	Sequence
D2-1 (SEQ ID NO: 13)	Forward	AGTTGTTAGTCTACGTGGACCGAC
D2-2389 (SEQ ID NO: 14)	Forward	GTGACACTGTATTGGGAGTCATGGTGCAG
D2-4996 (SEQ ID NO: 15)	Forward	GGTAATGGTGTGTTACAAGGAGTG
D2-5219 (SEQ ID NO: 16)	Forward	AAGCCCTTAGAGGACTTCCAATAA
CD2-373 (SEQ ID NO: 17)	Reverse	TCAGCATCCTTCCAATCTCTTTCC
CD2-2736 (SEQ ID NO: 18)	Reverse	CCGCAGAGATCGTTTTCCTGCCTG
cD2-5318 (SEQ ID NO: 19)	Reverse	CTAACTGGTGATAGCAGCCTCATGG
CD2-5358 (SEQ ID NO: 20)	Reverse	GAAATGGGCTTCGTCCATGATAATCAGG

TABLE 4

Primers for Sanger Sequencing Primer conc. = 10 uM in Tris buffer, pH 7.5 vRNA per Reaction = 100 ng in Tris buffer		
Primer	Reversion Region	Reversion nt site
D2-1 or cD2-373	5' NCR	57
D2-2389 or cD2-2736	NS1	2579
D2-4996 or cD2-5318	NS3	5270
D2-5219 or cD2-5358	NS3	5270

Pacbio Sequencing and Variant Analysis

RT-PCR to Amplify PCR Product

[0193] The optimized RT-PCR procedure was used to amplify cDNA from four TDV2 revertants viral RNA using gene specific primers (GSP). Reverse transcription was carried out on 6.5 µL of viral RNA with 3R reverse primer following by nested PCR using two best performing primer pairs (5F and 3R) at 25, 30, 35 and 40 PCR cycles. The optimum of 30 cycles was chosen and a 0.45× Ampure XP bead purification was used to remove non-specific amplification products less than 3 kilobases.

Pacbio Sequencing and Variant Analysis

[0194] The materials from RT-PCR reactions were used to construct PacBio sequencing libraries.

Materials and Methods for Pacbio's Amplicon Template Preparation and Sequencing Protocol:

Materials:

- [0195]** SMRTbell™ Template Prep Kit. Pacbio 100-259-100
[0196] DNA/Polymerase Binding Kit. Pacbio 100-372-700

- [0197]** MagBead Kit. Pacbio catalog number 100-676-500
[0198] DNA Sequencing Reagent. Pacbio 100-612-400
[0199] DNA Internal Control Complex. Pacbio 100-364-600
[0200] SMRT® Cells. Pacbio 100-171-800
[0201] Pre-washed AMPure XP Beads. Beckman Coulter A63881
[0202] Qiagen Buffer EB. Qiagen 19086
[0203] Magnetic rack. ThermoFisher Scientific 12321
[0204] DynaMag™-96 Bottom Magnet. ThermoFisher Scientific 12332D
[0205] Eppendorf™ DNA LoBind Microcentrifuge Tubes, Eppendorf 022431021
[0206] VWR® PCR 8-Well Tube Strips, VWR 20170-004
[0207] Qubit™ dsDNA BR Assay Kit, ThermoFisher Scientific Q32850
[0208] Agilent DNA 12000 Kit, Agilent 5067-1508
[0209] Hard-Shell PCR Plates 96-well, thin wall, Bio-Rad HSP9655

Pre-Wash AMPure XP Beads: (Use Washed AMPure XP Beads for all Library Prep Steps)

- [0210]** 1. Thoroughly re-suspend AMPure XP beads.
- [0211]** 2. Pipette 500 µl of AMPure XP beads into a 1.5 ml centrifuge tube.
- [0212]** 3. Spin for 1 minute at max speed in a bench top microcentrifuge.
- [0213]** 4. Place the tube on a magnetic rack until the supernatant clears.
- [0214]** 5. Transfer the supernatant to another 1.5 ml centrifuge tube and save for later.
- [0215]** 6. Add 1 ml of molecular biology grade water to the beads and vortex to re-suspend the beads completely.
- [0216]** 7. Centrifuge the beads for 1 minute at max speed in a bench top microcentrifuge.

- [0217] 8. Place the tube on a magnetic rack until the supernatant clears.
- [0218] 9. Discard the water.
- [0219] 10. Repeat the water wash 4 more times (steps 6-9).
- [0220] 11. Add 1 ml of Qiagen buffer EB to the bead pellet and vortex to completely re-suspend beads.
- [0221] 12. Centrifuge for 1 minute at max speed in a bench top microcentrifuge.
- [0222] 13. Place the bead tube on a magnetic rack until the supernatant clears.
- [0223] 14. Discard the EB buffer.
- [0224] 15. Completely re-suspend the beads in the original AMPure XP supernatant collected in step 5 by vortexing.
- [0225] 16. Store the washed beads at 4° C. for up to 3 months.

Purify PCR Product:

- [0226] 1. Determine the sample volume of your large-scale PCR reactions. In this case, 3×PCR reactions produced 150 µl of product for each sample. Move the PCR products to new 1.5 ml LoBind microcentrifuge tubes. Add 0.45× volumes (67.5 µl) of AMPure XP beads to each sample.
- [0227] 2. Mix the bead/DNA solution thoroughly.
- [0228] 3. Quickly spin down the tubes to collect the beads.
- [0229] 4. Allow the DNA to bind to the beads by rotating the tubes for 10 minutes at room temperature. Note that the bead/DNA mixing is critical to yield. After rotating, the bead/DNA mixture should appear homogenous. Failure to thoroughly mix the DNA with the bead reagent will result in inefficient DNA binding and reduced sample recoveries.
- [0230] 5. Centrifuge the tubes for 1-5 seconds.
- [0231] 6. Place the tubes in a magnetic bead rack until the beads collect to the side of the tubes and the solution appears clear. The actual time required to collect the beads to the side depends on the volume of beads added.
- [0232] 7. With the tubes still on the magnetic bead rack, slowly pipette off cleared supernatant. Avoid disturbing the bead pellet.
- [0233] 8. Wash the beads by slowly dispensing 1.5 ml of freshly prepared 70% ethanol against the side of the tube opposite the beads. Let the tubes sit for 30 seconds. Do not disturb the bead pellet. After 30 seconds, pipette and discard the ethanol.
- [0234] 9. Repeat step 8 above.
- [0235] 10. Remove residual 70% ethanol and dry the bead pellet. Remove tubes from the magnetic bead rack and spin to pellet beads. Both the beads and any residual 70% ethanol will be at the bottom of the tube. Place the tubes back on magnetic bead rack. Pipette off any remaining 70% ethanol.
- [0236] 11. Check for any remaining droplets in the tubes. If droplets are present, repeat step 10.
- [0237] 12. Remove the tubes from the magnetic bead rack and allow beads to air-dry with the tube caps open for 30 to 60 seconds.
- [0238] 13. Add 39 µl of Qiagen Buffer EB to the beads. Mix until homogenous by vortexing for 1 to 2 minutes at 2000 rpm. Spin the tubes down to pellet the beads

and then place the tubes on the magnetic rack until the supernatant clears. Collect 37 µl the eluted sample and place into new 200 µl PCR tubes. Place the sample on ice.

Repair DNA Damage:

- [0239] 1. Add the following reagents from the SMRTbell™ Template Prep Kit to each sample.

DNA Damage Repair Buffer	5 µl
NAD+	0.5 µl
ATP high	5 µl
dNTP	0.5 µl
DNA Damage Repair Mix	2 µl

Make a master mix of DNA Damage Repair Buffer, NAD+, ATP high, and dNTP when preparing more than 1 sample. Add 11 µl of master mix and 2 µl of DNA Damage Repair Mix to each sample.

- [0240] 2. Mix the reaction well by flicking the tube.
- [0241] 3. Centrifuge quickly to collect the reaction in the bottom of the tube.
- [0242] 4. Place the sample tubes in a thermocycler and incubate at 37° C. for 1 hr and then at 4° C. for 1 minute. Place the samples on ice.

Repair Ends:

- [0243] 1. Add the following reagents from the SMRTbell™ Template Prep Kit to each sample.

End Repair Mix	2.5 µl
----------------	--------

- [0244] 2. Mix the reaction well by flicking the sample tubes.
- [0245] 3. Centrifuge briefly to collect samples.
- [0246] 4. Place the samples in a thermocycler and incubate at 25° C. for 5 minutes and then place the samples back on ice. Proceed immediately to the next step and do not allow the reaction to incubate longer than 5 minutes.

Purify End Repaired Product:

- [0247] 1. Add 0.45× volumes (23.6 µl) of AMPure XP beads to each sample.
- [0248] 2. Mix the bead/DNA solution thoroughly by flicking the tubes.
- [0249] 3. Quickly spin down the tubes to collect the beads.
- [0250] 4. Allow the DNA to bind to the beads by rotating the tubes for 10 minutes at room temperature. Note that the bead/DNA mixing is critical to yield. After rotating, the bead/DNA mixture should appear homogenous. Failure to thoroughly mix the DNA with the bead reagent will result in inefficient DNA binding and reduced sample recoveries.
- [0251] 5. Centrifuge the tubes for 1-5 seconds.
- [0252] 6. Place the tubes in a DynaMag™-96 Bottom Magnet bead rack until the beads collect to the side of the tube and the solution appears clear.

- [0253] 7. With the tubes still on the magnetic bead rack, slowly pipette off cleared supernatant. Avoid disturbing the bead pellet.
- [0254] 8. Wash the beads by slowly dispensing 200 μ l of freshly prepared 70% ethanol against the side of the tube opposite the beads. Let the tubes sit for 30 seconds. Do not disturb the bead pellet. After 30 seconds, pipette and discard the ethanol.
- [0255] 9. Repeat step 8 above.
- [0256] 10. Remove residual 70% ethanol and dry the bead pellet. Remove tubes from the magnetic bead rack and spin to pellet beads. Both the beads and any residual 70% ethanol will be at the bottom of the tubes. Place the tubes back on magnetic bead rack. Pipette off any remaining 70% ethanol.
- [0257] 11. Check for any remaining droplets in the tubes. If droplets are present, repeat step 10.
- [0258] 12. Remove the tubes from the magnetic bead rack and allow beads to air-dry with the tube caps open for 30 to 60 seconds.
- [0259] 13. Add 33 μ l of Qiagen Buffer EB to the beads. Mix until homogenous by vortexing for 1 to 2 minutes at 2000 rpm. Spin the tubes down to pellet the beads and then place the tubes on the magnetic rack until the supernatant clears. Collect 31 μ l of each eluted sample and place into new 200 μ l PCR tubes.
- [0260] 14. Place the samples on ice.

Blunt Ligation:

- [0261] 1. Add the following reagents from the SMRTbell™ Template Prep Kit to each sample. If making a master mix, combine all of the components below except for the Blunt Adapter. Aliquot 7 μ l of the master mix to each sample and then add the adapter directly to the mixture.

Blunt Adapter	2 μ l
Template Prep Buffer	4 μ l
ATP Lo	2 μ l
Ligase	1 μ l

- [0262] 2. Mix the reaction by flicking the tubes.
- [0263] 3. Centrifuge the tubes briefly to collect the samples.
- [0264] 4. Place the samples in a thermocycler and incubate at 25° C. for 1 hr and then at 65° C. for 10 minutes to inactivate the ligase. Place the samples back on ice.
- [0265] 5. Add 0.5 μ l of Exo III and 0.5 μ l ExoVII to each sample and incubate the samples at 37° C. for 1 hr.
- [0266] 6. You must proceed to the next step.

Purify SMRTbell Libraries

- [0267] 1. Add 0.45 $\times$ volumes (18.45 μ l) of AMPure XP beads to each sample.
- [0268] 2. Mix the bead/DNA solution thoroughly by flicking the tubes.
- [0269] 3. Quickly spin down the tubes to collect the beads.
- [0270] 4. Allow the DNA to bind to the beads by rotating the tubes for 10 minutes at room temperature. Note that the bead/DNA mixing is critical to yield. After rotating, the bead/DNA mixture should appear

homogenous. Failure to thoroughly mix the DNA with the bead reagent will result in inefficient DNA binding and reduced sample recoveries.

- [0271] 5. Centrifuge the tubes for 1-5 seconds.
- [0272] 6. Place the tubes in DynaMag™-96 Bottom Magnet bead rack until the beads collect to the side of the tube and the solution appears clear.
- [0273] 7. With the tubes still on the magnetic bead rack, slowly pipette off cleared supernatant. Avoid disturbing the bead pellet.
- [0274] 8. Wash the beads by slowly dispensing 200 μ l of freshly prepared 70% ethanol against the side of the tube opposite the beads. Let the tubes sit for 30 seconds. Do not disturb the bead pellet. After 30 seconds, pipette and discard the ethanol.
- [0275] 9. Repeat step 8 above.
- [0276] 10. Remove residual 70% ethanol and dry the bead pellet. Remove tubes from the magnetic bead rack and spin to pellet beads. Both the beads and any residual 70% ethanol will be at the bottom of the tubes. Place the tubes back on the magnetic bead rack. Pipette off any remaining 70% ethanol.
- [0277] 11. Check for any remaining droplets in the tubes. If droplets are present, repeat step 10.
- [0278] 12. Remove the tubes from the magnetic bead rack and allow beads to air-dry with the tube caps open for 30 to 60 seconds.
- [0279] 13. Add 52 μ l of Qiagen Buffer EB to the beads. Mix until homogenous by vortexing for 1 to 2 minutes at 2000 rpm. Spin the tubes down to pellet the beads and then place the tubes on the magnetic rack until the supernatant clears. Collect 50 μ l of each eluted sample and place into new 200 μ l PCR tubes.
- [0280] 14. Add 0.45 $\times$ volumes (22.5 μ l) of AMPure Xp beads to the elute sample. Repeat steps 2-12.
- [0281] 15. Add 12 μ l of Qiagen Buffer EB to the beads. Mix until homogenous by vortexing for 1 to 2 minutes at 2000 rpm. Spin the tubes down to pellet the beads and then place the tubes on the magnetic rack until the supernatant clears. Collect 10 μ l of the eluted library.
- [0282] 16. Quantify 1 μ l of sample library using a Qubit™ dsDNA BR Assay Kit and the Qubit 3.0 Fluorometer
- [0283] 17. Determine the size distribution of the library by running 1 μ l (0.5-50 ng/ μ l) of sample library on a 2100 Bioanalyzer using an Agilent DNA 12000 Kit.

Anneal and Bind SMRTbell Libraries

[0284] The DNA/Polymerase Binding and Magbead kits are required for this step. Use the Pacbio P6C4 Binding Calculator, version 2.3.1.1, to determine the conditions for annealing the sequencing primer and binding polymerase to SMRTbell libraries. See the Binding Calculator parameters below.

Sample Name	P1	P3	P5	P51
Sample Volume to Use	3 μ l	3 μ l	3 μ l	3 μ l
# of SMRTCells				
Concentration	23.9 ng/ μ l	15.7 ng/ μ l	39.4 ng/ μ l	26.9 ng/ μ l
Insert size	6300 bp	6300 bp	6300 bp	6300 bp

-continued

Sample Name	P1	P3	P5	P51
Size Selection Protocol	No MagBead OCPW	No MagBead OCPW	No MagBead OCPW	No MagBead OCPW
Binding Kit	P6	P6	P6	P6
Preparation	Small	Small	Small	Small
Long-Term Storage	No	No	No	No
DNA Control	No	No	No	No
Complex Reuse	No	No	No	No
Standard	Yes	Yes	Yes	Yes
Concentration	Custom	Custom	Custom	Custom
OnPlate	(0.0125 nM)	(0.0125 nM)	(0.0125 nM)	(0.0125 nM)
Control To	Default (1.2%)	Default (1.2%)	Default (1.2%)	Default (1.2%)
Template	Default (10)	Default (10)	Default (10)	Default (10)
Polymerase: Template Ratio	Default (20)	Default (20)	Default (20)	Default (20)

[0285] 1. Dilute the Sequencing Primer by adding the reagents below into a new 200 P3 PCR tube.

[0286] a. Sequencing Primer 1 μ l

[0287] b. Elution Buffer 32.3 μ l

[0288] 2. Incubate at 80° C. for 2 minutes then hold at 4° C. Conditioned sequencing primer may be stored at -20° C. and used for up to 30 days.

[0289] 3. Annealing the diluted sequencing primer by adding the appropriate amount of reagents in the order listed below into 200 μ l PCR tubes.

Sample Name	P1	P3	P5	P51
Volume H2O	13.6 μ l	7.9 μ l	24.3 μ l	15.7 μ l
10x Primer Buffer	2.1 μ l	1.4 μ l	3.5 μ l	2.4 μ l
Sample Volume	3 μ l	3 μ l	3 μ l	3 μ l
Diluted Sequencing Primer	2.3 μ l	1.5 μ l	3.8 μ l	2.6 μ l
Total Volume	21 μ l (at 0.8333 nM)	13.8 μ l (at 0.8333 nM)	34.6 μ l (at 0.8333 nM)	23.6 μ l (at 0.8333 nM)

[0290] 4. Incubate at 20° C. for 30 minutes.

[0291] 5. Place the samples on ice.

[0292] 6. Dilute SA-DNA Polymerase by adding the reagents below to a new 200 μ l PCR tube.

a.	SA-DNA Polymerase P6	3 μ l
b.	Binding Buffer v2	27 μ l
c.	Total Volume (50 nM final)	30 μ l

[0293] 7. Add the following components to the primer-annealed libraries and mix by flicking the tubes. Spin down to collect the sample.

Sample Name	P1	P3	P5	P51
dNTP	3.5 μ l	2.3 μ l	5.8 μ l	3.9 μ l
DTT	3.5 μ l	2.3 μ l	5.8 μ l	3.9 μ l

-continued

Sample Name	P1	P3	P5	P51
Binding Buffer v2	3.5 μ l	2.3 μ l	5.8 μ l	3.9 μ l
Polymerase Dilution	3.5 μ l	2.3 μ l	5.8 μ l	3.9 μ l
Total Volume	35 μ l (at 0.5 nM)	23 μ l (at 0.5 nM)	57.7 μ l (at 0.5 nM)	39.4 μ l (at 0.5 nM)

[0294] 8. Close the PCR tubes and incubate at 30° C. for 30 minutes and then return to ice. The polymerase bound SMRTBell libraries can be stored at 4° C. for ~7 days and should be sequenced within 10 days. Sequencing yield and quality begins to decline after 7 days.

[0295] 9. Prepare the SMRTBell libraries for sequencing by mixing the following for each sample in new 200 μ l PCR tubes.

a.	MagBead Binding Buffer	8.8 μ l
b.	Polymerase Bound Library	0.22 μ l

[0296] 10. Place the diluted libraries back on ice.

[0297] 11. Mix the MagBeads by vortexing until they are homogenous.

[0298] 12. For each sample, add 35 μ l of Magbeads to an empty 1.5 ml LoBind tube. For more than one sample make a master mix and add 10% to account for pipetting loss. In this case, we have four samples so $4.1 \times 35 \mu\text{l} = 143.5 \mu\text{l}$.

[0299] 13. Centrifuge for a couple of seconds to collect the MagBeads and place them on a magnetic bead rack until the supernatant clears.

[0300] 14. Discard the supernatant.

[0301] 15. Add 143.5 μ l of MagBead Wash Buffer to the beads and mix until homogenous by pipetting 10 times.

[0302] 16. Centrifuge the beads for a couple of seconds to collect the beads and place back on the magnetic rack until the supernatant clears.

[0303] 17. Discard the supernatant.

[0304] 18. Add 143.5 μ l of MagBead Binding Buffer to the beads and mix by pipetting 10 times.

[0305] 19. Centrifuge the beads for a couple of seconds to collect the beads.

[0306] 20. Aliquot 35 μ l of beads to four separate 1.5 ml LoBind tubes and place on ice.

[0307] 21. Place the beads on a magnetic bead rack until the supernatant clears and then discard the supernatant.

[0308] 22. Add 9 μ l of the diluted library, from step 10, to the corresponding bead tube. Cap the bead tubes and mix the beads until homogenous by flicking the tubes. Centrifuge for a couple of seconds to collect the beads and then place the beads back on ice.

[0309] 23. Incubate the beads at 4° C. on a rotator for 30 minutes and then place back on ice.

[0310] 24. While the beads are incubating, dilute the DNA Internal Control Complex in new 200 μ l PCR tubes. Mix the following.

[0311] a. Control 1 Mix

i.	MagBead Binding Buffer	99 μ l
ii.	DNA Internal Control Complex	1 μ l

[0312] b. Control 2 Mix

i.	MagBead Binding Buffer	49 μ l
ii.	Control 1 mix	1 μ l

[0313] 25. Centrifuge the beads from step 23 for a couple of seconds and then place on a magnetic bead rack until the supernatant clears. Discard the supernatant and place the beads back on ice.

[0314] 26. Add 18 μ l of MagBead Binding buffer to the beads and mix by flicking the tubes.

[0315] 27. Place the beads on a magnetic bead rack until the supernatant clears and then discard the supernatant.

[0316] 28. Add 18 μ l of MagBead Wash buffer to the beads and mix by flicking the tubes.

[0317] 29. Place the beads on a magnetic bead rack until the supernatant clears and then discard the supernatant.

[0318] 30. Add 3 μ l of the diluted Control 2 Mix and 42 μ l of MagBead Binding buffer to the beads and mix by flicking the tubes. Centrifuge for a couple of seconds to collect the bead mixture and place on ice until use.

[0319] 31. Place a new 96-well PCR plate on ice.

[0320] 32. Mix the beads from step 30 until homogenous by pipetting and aliquot 45 μ l of each of the four mixtures to separate wells in the same column (e.g., column 1, rows A-D) of the 96-well PCR plate that is on ice.

Pacbio Loading and Data Collection Time

[0321] 1. The Pacbio run was configured using Pacbio RS Remote software, version 2.3. The following settings were used:

[0322] a. Collection Protocol=MagBead OneCellPerWell v1

[0323] b. AcquisitionTime (Movie Time)=360

[0324] c. InsertSize=5000

[0325] d. StageHS=True

[0326] e. SizeSelectionEnabled=False

[0327] f. DNAControlComplex=2 kb_control

[0328] g. Use2ndLook=False

[0329] h. NumberOfCollections=1

[0330] 2. The PCR Plate containing the sequencing libraries, DNA Sequencing Reagents and SMRTCells were loaded into the Pacbio RS II. The libraries were sequenced for 360 minutes.

Variant Detection and Quantification

[0331] The raw sequencing data and dengue virus reference genome (GenBank accession U87411; Dengue Virus type 2, Strain 16681) were imported into SMRTanalysis software version 2.3.0.140936.p5.167094 that consisted of Daemon 2.3.0.139497, SMRTpipe 2.3.0.139497, SMRT Portal 2.3.0.140893 and SMRT View 2.3.0.140836. The P1, P3, P5 and P51 reads were aligned to the dengue virus reference genome (GenBank accession U87411; Dengue Virus type 2, Strain 16681) using the RS_ReadsOfinsert_Mapping.1 protocol of the Smrtanalysis software with the default mapping parameters and the following filtering parameters, Minimum Full Passes=5, Minimum Predicted Accuracy=90, Minimum Read Length of Insert=5000, Maximum Read Length of Insert=7000. The reads of insert fasta (fasta and fastq files can be used interchangeably) files generated by the RS_ReadsOfinsert_Mapping.1 protocol were aligned to the dengue virus (GenBank accession U87411; Dengue Virus type 2, Strain 16681) reference genome with BWA and variants at each attenuation site were identified from the alignment SAM files using the Sam2Tsv (<http://lindenb.github.io/jvarkit/Sam2Tsv.html>) and custom bash scripts. A custom perl script was used to quantify variants at the three attenuation sites of interest.

[0332] Four SMRTcells were sequenced, producing average about 8.7 kb reads of insert post filtration. The reads were aligned to the wildtype dengue virus reference genome (GenBank accession U87411) using BLASR and variants were called with QUIVER. The data are summarized in the Result section.

Further RT-PCR Optimization to Produce a Single PCR Product (Free of Non-Specific Bands) Using TDV2

[0333] The existing protocol from phase I did not produce a single and clean amplicon PCR product. Therefore, additional protocol optimizations and methods were attempted to produce a PCR product with a single amplicon. These optimizations included using a different combination of previously designed primers for reverse transcription and PCR, increasing the RT and touchdown PCR reaction temperatures. These various optimization steps are described below:

Reverse Transcription and Touchdown PCR Optimization

[0334] To identify new combinations PCR primers for optimization, in silico or ePCR with previously designed primer sets was done, allowing for 30% primer binding site mismatch, with Unipro UGENE software. Based on the ePCR analysis that demonstrated low off target binding and our results with nested PCR, primers 5F and 3R were selected for optimization. Primer sequences are in Table 5.

Two Pairs of Primers Used for RT-PCR Optimization in Nested PCR Reaction

ber of UMIs in the reads rather than total reads, since total reads are skewed due to non-uniform amplification. After

TABLE 5

Primer sequences and characteristics. Note that the 1R reverse primer is also the prime used for reverse transcription. The primers are in pairs.						
Name	Dir.	Sequence	Start	Stop	Tm	Product
1F (SEQ ID	F	AAGACAGATTCTTTGAGGGAGCTA	26	49	59	5368
1R (SEQ ID	R	GGTCTGTGAAATGGGCTTCG	5393	5374	59	5368
5F (SEQ ID	F	GACAGATTCTTTGAGGGAGCTAAG	28	51	59	5364
3R (SEQ ID	R	GTCTGTGAAATGGGCTTCGTC	5392	5372	60	5364

[0335] First, reverse transcription was done using 6.5 µl of TDV2 viral RNA (~10⁸ viral genome copies/µl) with primer 3R at two different reaction temperatures 50° C. and 55° C. Then, touchdown PCR was done with primers 5F and 3R as follows: 95° C. for 30 sec, annealing at 68° C. for 20 sec with a 1 degree drop each cycle for 10 cycles (i.e., cycle 1-68° C., cycle 2-67° C., cycle 3-65° C. . . .) and extension at 68° C. for 6 min. An additional 5 cycles were run at 95° C. for 30 sec, annealing at 53, 55.4, 57.2, 59.4, 61.9, 65.7, or 68° C. for 20 sec and extension at 68° C. for 6 min, giving a total of 15 cycles. PCR product was run in a 1% TAE agarose gel and the results show that a single PCR amplicon of expected size was produced.

primers and RT-PCR conditions were optimized, the next step was to incorporate UMIs into the 3R reverse transcription primer and thus, into each viral molecule to improve the quantification of the observed variants.

[0337] A fully random 16-nucleotide UMI and a 5' PCR Primer IIA binding site were added to the 5' end of the 3R primer. Primer sequences and characteristics are shown in Table 6.

Primers and Oligonucleotide Sequences Used During Incorporation of UMIs

TABLE 6

Primer and oligonucleotide sequences and characteristics. 5F and UMI 3R were used to introduce UMI using the previously established successful amplification using 5F and 3R. SMARTER IIA and 5' PCR Primer IIA were used with UMI 3R when the template switching protocol was attempted.						
Name	Sequence	Start	Stop	Tm	Product	
5F	GACAGATTCTTTGAGGGAGCTAAG	28	51	59	5405	
5' PCR Primer IIA	AAGCAGTGGTATCAACGCAGAGTAC	28	51	59	5405	
SMARTER IIA	AAGCAGTGGTATCAACGCAGAGTAC XXXXX				5462	
UMI 3R	AAGCAGTGGTATCAACGCAGAGTAC NNNNNNNNNNNNNNNGTCTGTGAA ATGGGCTTCGTC					

X = undisclosed nucleotide based on proprietary sequence from Clontech Laboratories, Inc.

Introduction of Unique Molecular Identifiers

[0336] Unique molecular Identifiers (UMIs) are advantageous in that they can correct for artifacts produced by library preparation, such as PCR duplication and preferential amplification of certain variants. Also, UMIs have been shown to improve detection of rare variants by correcting for sequencing errors. The concept of UMI is that each RNA or DNA molecule is attached with a molecular tag or barcode. Amplification products showing same UMIs are duplication of same DNA molecule. Sequence reads having different UMIs are the original molecules. UMIs have the potential to increase the accuracy of detection of mutations at 1% frequency or lower by removing false positives. The target quantification can be better achieved by counting the num-

TABLE 7

Primers and SEQ ID numbering	
Name	Direction
5F (SEQ ID NO: 9)	F
5' PCR Primer IIA (SEQ ID NO:21)	Both
SMARTER IIA	5' RT Adapter
UMI 3R (SEQ ID NO: 22)	3' RT Adapter containing UMI

[0338] Next, RT-PCR was performed with the optimized conditions using 6.5 µl of TDV2 viral RNA (~10⁸ viral genome copies/µl). The PCR product was run in a 1% TAE

agarose gel. Subsequently, large scale PCR was performed. A 0.45× Ampure XP bead purification was used to concentrate the PCR product and remove excess primers. Large scale PCR and Ampure XP bead purification gave strong amplification of a single band at the expected size.

Pacbio Sequencing and Variant Analysis

[0339] The material from optimized RT-PCR reaction was used to construct a PacBio sequencing library using Pacbio's amplicon library prep protocol as described above.

Quantification of TDV2 Genomes Using Read Based and UMI Based Methods

[0340] Variants were not called with QUIVER as they were in Phase I. The reason is that typical variant calling algorithms do not maintain individual read ID information since variants are called based on read consensus and coverage at variant sites. Therefore, individual read information is not needed. However, to utilize UMI based quantification methods, multiple reads containing the same UMI and sequence are collapsed into a single instance of that read. Consequently, it is necessary to maintain individual read ID information when using UMI based quantification methods and as such, typical variant callers do not facilitate analysis with UMIs. Therefore, an alternative variant identification workflow was developed to facilitate the use of UMIs. Briefly, this workflow consisted of trimming off the SMARTER UMI 3R adapter sequence, extracting the UMI sequence and appending it to read ID of the TDV2 read of insert fastq file, aligning the trimmed fastq file to the Dengue virus reference genome (GenBank accession U87411; Dengue Virus type 2, Strain 16681) using BWA-mem and counting the revertants at each attenuation site.

Apply Optimized Library Prep and Variant Detection Protocol to TDV2 Mono-Revertants

Large Scale PCR to Incorporate UMIs into Mono-Revertants

[0341] Large scale PCR was performed as described above in the "Introduce Unique Molecular Identifiers into Each Viral Molecule" section. The PCR products were run in a 1% TAE agarose gel. The results show a single amplicon of expected size for each revertant.

Pacbio Sequencing and Variant Analysis

[0342] The material was used to construct Pacbio sequencing libraries using Pacbio's amplicon library prep protocol as outlined above. One SMRTcell was sequenced per revertant.

[0343] Next, variant analysis was done as described in the "Quantify TDV2 Variants Using Read Based And UMI Based Methods" section.

Apply Optimized Library Prep and Variant Detection Protocol to Four Complex Revertant Mixtures

Preparation of Complex Revertant Mixtures

[0344] Four complex revertant mixture pools were prepared based on the genome copies in each individual revertant as determined by qRT-PCR and transported to the

analytical group without disclosing any information on the composition. Details of the pools are described in the result section.

Large Scale PCR to Incorporate UMIs into Complex Revertant Mixtures

[0345] RT was performed on four complex revertant mixtures, Complex Mix 1-4, as described above in the "Introduce Unique Molecular Identifiers into Each Viral Molecule" section using 6.5 µl of viral RNA (~2.8×10⁷ viral genome copies/µl) and the SMARTER UMI 3R primer. Next, large scale touchdown PCR was performed as described above and PCR products were run in a 1% TAE agarose gel. The results show a single amplicon of expected size for Complex Mix 1-4.

Pacbio Sequencing and Variant Analysis

[0346] The material was used to construct PacBio sequencing libraries using Pacbio's amplicon library prep protocol as described above. One SMRTcell was sequenced per complex revertant mixture.

[0347] Next, variant analysis was done as described in the "Quantify TDV2 Variants Using Read Based and UMI Based Methods" section.

Results

Optimization of Pacbio Sequencing for TDV2 Analysis

RT-PCR Optimization

[0348] In touchdown PCR reactions, there was a band present in all lanes in the expected size range, plus several non-specific amplification products.

[0349] Two best performing primer pairs (1 and 4) optimized for number of cycles in nested PCR gave strong amplification of an expected size band. Primer pair 1 had fewer non-specific products than pair 4 and was taken forward to annealing temperature optimization. The optimum of 11 cycles was chosen and a 0.45× Ampure XP bead purification was used to remove non-specific amplification products less than 3 kilobases.

PacBio Sequencing and Variant Analysis

[0350] The material was used to construct a PacBio sequencing library as described above. One SMRTcell was sequenced, yielding 366 Mbp of data and 12,467 reads of insert. The Mean number of passes around each read of insert was 6.15 giving a mean read quality of 0.9828. The mean read of insert length was 5321 nucleotides.

[0351] Note: The mean number of passes is the mean number of times a polymerase read the same molecule. Derived by the PacBio analytics software, PacBio sequencing templates are single stranded circles. The polymerase keeps going around until it falls off becomes inactivated by the laser, or the run finishes.

[0352] The alignment of reads to the dengue virus reference genome (GenBank accession U87411; Dengue Virus type 2, Strain 16681) initially identified eleven variants. The mapped reads of insert length are shown in FIG. 4.

[0353] The first three of these were poorly supported by low coverage were discarded from the analysis. The remaining variants were well supported by the percentage of reads calling each one. The data are summarized in Table 8.

TABLE 8

Summary of the observed variants with the numbers of supporting reads.				
Position	Variant	Reads total	Reads calling	Percentage
57	C > T	11,917	11,846	99.4
524	A > T	11,121	11,080	99.6
592	A > G	12,236	12,212	99.8
900	T > C	12,403	11,853	95.6
2,055	C > T	11,339	11,297	99.6
2,579	G > A	11,710	11,461	97.9
4,018	C > T	11,849	11,812	99.7
5,270	A > T	11,983	11,882	99.1

Table 8:

[0354] The combinations of the three attenuation variants at positions, 57, 2579 and 5270 were examined and are summarized in Table 9.

TABLE 9

Haplotype analysis of the three attenuation loci.				
Count	Percent	57	2579	5270
10472	96.703	T	A	T
28	0.259	C	A	T
15	0.139	A	A	T
12	0.111	G	A	T
201	1.856	T	G	T
5	0.046	T	C	T
48	0.443	T	a	A
21	0.194	T	A	C
17	0.157	T	A	G
2	0.018	C	G	T
1	0.009	A	G	T
4	0.037	T	G	A
1	0.009	T	G	G
1	0.009	C	A	C
1	0.009	A	A	C
0	0	C	G	A

[0355] The read counts are shown in the first column. The data are arranged such that observed variants are in order of positions 57, 2579 and 5270. Although PCR limits the ability to infer quantitative information about the input RNA, the TDV2 haplotype TAT is by far the most abundant. No reads having the reference haplotype CGA were found. The remainder have either one or two variants with respect to the TDV2 haplotype.

Identification of Mono Revertants from Blinded Samples

Sanger Sequencing to Identify Mono Revertants

[0356] Following conditions were selected as final optimized conditions:

[0357] RT at 55° C. for 30 minutes

[0358] Annealing Temp.: 60° C.

[0359] Number of PCR cycles: 30

[0360] Gel extraction kit from Qiagen was used to purify all PCR products of all TDV2 revertants. Sanger sequencing analysis confirmed the identity of TDV2 revertants (Table 10).

TABLE 10

Summary of Sanger Sequencing Results		
Revertants	Reversion Site	Confirmed (Yes/No)
P1	A—NS1-2579—G	Yes
P3	T—NS3-6270—A	Yes
P5	T—5'NC-57—C	Yes
P51	T—5'NC-57—C	Yes
	A—NS1-2579—G	
TDV2	None-No Reversion	Yes

Pac Bio Sequencing to Identify Mono Revertants

Analysis of Mono Revertants

[0361] 30 cycles were determined to be optimal because the PCR product began to smear up the gel at 35-40 cycles, which indicated PCR chimera.

Pacbio Sequencing and Variant Analysis

[0362] The Pacbio sequencing statistics are shown in Table 11. The P1, P3, P5 and P51 reads were aligned to the dengue virus reference genome (GenBank accession U87411; Dengue Virus type 2, Strain 16681) using BLASR. The numbers of mapped reads are shown in Table 11. Variants were called with QUIVER. Haplotypes T:57-G:2579-T:5270, T:57-A:2579-A:5270, C:57-A:2579-T:5270 and C:57-G:2579-T:5270 were the most abundant haplotypes, making up >99% of the population, in P1, P3, P5 and P51, respectively (Table 12).

TABLE 11

Pacbio sequencing and mapping statistics for mono-revertants, P1, P3, P5 and P51.						
Sample	Yield (Mbps)	Reads of Insert	Mean number of Passes	Mean Read Quality	Mean Read of Insert Length (bps)	Mapped Reads
P1	49	9232	8.89	0.9955	5343	9240
P3	29	5590	8.97	0.9955	5341	5596
P5	59.7	11182	8.38	0.9954	5344	11199
P51	47.9	8965	8.65	0.9957	5347	8975

TABLE 12

Variant counts at three attenuation sites.		
P1 ATTENUATION SITE	READ COUNT	READ PERCENT
T:57-G:2579-T:5270	8476	99.42522
T:57-G:2579-A:5270	17	0.199413
T:57-A:2579-T:5270	12	0.140762
T:57-G:2579-C:5270	11	0.129032
C:57-G:2579-T:5270	9	0.105572
P2 ATTENUATION SITE	READ COUNT	READ PERCENT
T:57-A:2579-A:5270	5232	99.24127
T:57-G:2579-A:5270	40	0.758725

TABLE 12-continued

P3 ATTENUATION SITE	READ COUNT	READ PERCENT
②	②	②
P51	READ	READ
ATTENUATION SITE	COUNT	PERCENT
C:57-G:2579-T:5270	8259	99.24297
C:57-A:2579-T:5270	28	0.336458
C:57-G:2579-C:5270	19	0.228311
C:57-G:2579-A:5270	16	0.192261

② indicates text missing or illegible when filed

TABLE 13

Summary of revertant analysis using conventional RT-PCR and pacbio sequencing					
Clone-derived viruses (Revertants)	Reversion site	5'NCR (nt57)	NS1 (nt2579)	N53 (nt5270)	
D2/IC-P1	NS1	C = 0.10% WT	G = 99.86% WT, A = 0.14%	Vac A = 0.20%	WT
D2/IC-P3	NS3	C = 0.00% WT	G = 0.75% WT	A = 100%	WT
D2/IC-P5	5'NCR	C = 100% WT	G = 0.57% WT	A = 0.17%	WT
D2/IC-P51	5'NCR & NS1	C = 100% WT	G = 99.66% WT, A = 0.34%	Vac A = 0.19%	WT
TDV2 (Vaccine) Raverants					
D2 = DENV2 16681		5'NCR-nt57	T	C	
IC = Infectious Clone	Key	NS1-nt2579	A	G	
Px = Parental (where x = 5, 1, and/or 3)		NS3-nt6270	T	A	

Further RT-PCR Optimization to Produce a Single PCR Product (Using TDV2 Viral RNA)

Reverse Transcription and Touchdown PCR Optimization

[0363] Further optimization of RT and touchdown PCR reaction produced a single PCR amplicon of expected size. The results show that a single PCR amplicon of expected size was produced. Based on these results, 55° C. was determined to be the optimal RT temperature since there was less banding and smearing below 5 kbps. 57° C. was determined to be the optimal annealing temperature for the last five cycles of the PCR since the amplicon abundance begins to decline above 57.2° C.

Introduction of Unique Molecular Identifiers

[0364] The UMI sequences were added to gene specific reverse 3R primer. After adding UMIs, the RT-PCR was performed with the optimized conditions using 6.5 µl of TDV2 viral RNA (~10⁸ viral genome copies/µl). The PCR product was run in a 1% TAE agarose gel. Subsequently, large scale PCR was performed. A 0.45× Ampure XP bead purification was used to concentrate the PCR product and remove excess primers. Large scale PCR and Ampure XP bead purification gave strong amplification of a single band.

Pacbio Sequencing and Variant Analysis

[0365] One SMRTcell was sequenced and yielded 42.9 Mbps of data and 8038 reads of insert. The mean number of

passes around each read of insert was 4.99 giving a mean read quality of 0.9838. The mean read of insert length was 5346 nucleotides.

Quantify TDV2 Genomes Using Read Based and UMI Based Methods

[0366] The variant counts using reads and UMIs are shown in Table 14. The results indicate that the T:57-A:2579-T:5270 TDV2 haplotype is by far the most common, while the wild-type haplotype (C:57-G:2579-A:5270) was not observed. However, some revertant haplotypes were present, including P1 (T:57-G:2579-T:5270), P3 (T:57-A:2579-A:5270) and P5 (C:57-A:2579-T:5270).

TABLE 14

TDV2 variant counts using reads or UMIs at three attenuation sites.				
ATTENUATION SITE	READ COUNT	READ PERCENT	UMI COUNT	UMI PERCENT
T:57-A:2579-T:5270	7041	96.17538588	6698	96.15274189
T:57-G:2579-T:5270	184	2.513317853	177	2.540913006
T:57-A:2579-G:5270	33	0.450758093	31	0.445018662
G:57-A:2579-T:5270	18	0.245868051	17	0.244042492
C:57-A:2579-T:5270	14	0.191230706	13	0.186620729
A:57-A:2579-T:5270	9	0.122934025	8	0.114843526
T:57-A:2579-C:5270	8	0.109274689	8	0.114843526
T:57-C:2579-T:5270	6	0.081956017	6	0.086132644
T:57-A:2579-A:5270	3	0.040978008	3	0.043066322

Apply Optimized Library Prep and Variant Detection Protocol to TDV2 Mono-Revertants

Large Scale PCR to Incorporate UMIs into Mono-Revertants

[0367] At this stage, the RT-PCR and library prep protocols were optimized. Therefore, we skipped the PCR optimization steps and proceeded directly to large scale PCR as described above in the “Introduce Unique Molecular Identifiers into Each Viral Molecule” section. The PCR products were run in a 1% TAE agarose gel. The results show a single amplicon of expected size for each revertant.

Pacbio Sequencing and Variant Analysis

[0368] One SMRTcell was sequenced per revertant. Next, variant analysis was done as described in the “Quantify TDV2 Variants Using Read Based and UMI Based Methods” section. The results are shown in the following tables of Table 15 (1 of 4 to 4 of 4).

Table 15: Variant Counts Tables Using Reads or UMIs at Three Attenuation Loci for TDV2 Mono Revertants (P1, P3, P5 and P51)

TABLE 15

(1 of 4): P1 variant counts using reads or UMIs at three attenuation sites.				
ATTENUATION SITE	READ COUNT	READ PERCENT	UMI COUNT	UMI PERCENT
T:57-G:2579-T:5270	11993	99.38675727	11121	99.39226026
C:57-G:2579-T:5270	30	0.248611917	29	0.259183126
T:57-G:2579-C:5270	13	0.107731831	13	0.116185539
T:57-A:2579-T:5270	12	0.099444767	11	0.098310841
T:57-G:2579-G:5270	7	0.058009447	7	0.062561444
G:57-G:2579-T:5270	7	0.058009447	5	0.044686746
T:57-T:2579-T:5270	2	0.016574128	2	0.017874698
T:57-G:2579-A:5270	2	0.016574128	2	0.017874698
A:57-G:2579-T:5270	1	0.008287064	-1	-0.008937349

TABLE 15

(2 of 4): P3 variant counts using reads or UMIs at three attenuation sites.				
ATTENUATION SITE	READ COUNT	READ PERCENT	UMI COUNT	UMI PERCENT
T:57-A:2579-A:5270	7924	98.98813242	7542	98.97637795
T:57-G:2579-A:5270	52	0.649594004	49	0.643044619
C:57-A:2579-A:5270	11	0.137414116	11	0.144356955
T:57-A:2579-G:5270	4	0.04996877	4	0.052493438
T:57-T:2579-A:5270	4	0.04996877	4	0.052493438
T:57-C:2579-A:5270	3	0.037476577	3	0.039370079
T:57-A:2579-C:5270	2	0.024984385	2	0.026246719
T:57-A:2579-T:5270	2	0.024984385	2	0.026246719
A:57-G:2579-A:5270	1	0.012492192	1	0.01312336
G:57-A:2579-A:5270	1	0.012492192	1	0.01312336
T:57-G:2579-C:5270	1	0.012492192	1	0.01312336

TABLE 15

(3 of 4): P5 variant counts using reads or UMIs at three attenuation sites.				
ATTENUATION SITE	READ COUNT	READ PERCENT	UMI COUNT	UMI PERCENT
C:57-A:2579-T:5270	9759	98.58571573	9177	98.6349957
C:57-G:2579-T:5270	79	0.79806041	72	0.773860705
C:57-A:2579-C:5270	33	0.333367007	29	0.311693895
C:57-A:2579-G:5270	10	0.101020305	9	0.096732588
T:57-A:2579-T:5270	7	0.070714214	6	0.064488392
A:57-A:2579-T:5270	3	0.030306092	3	0.032244196
C:57-A:2579-A:5270	3	0.030306092	3	0.032244196
C:57-C:2579-T:5270	2	0.020204061	2	0.021496131
C:57-T:2579-T:5270	1	0.010102031	1	0.010748065
G:57-A:2579-T:5270	1	0.010102031	1	0.010748065
G:57-G:2579-T:5270	1	0.010102031	1	0.010748065

TABLE 15

(4 of 4): P51 variant counts using reads or UMIs at three attenuation sites.				
ATTENUATION SITE	READ COUNT	READ PERCENT	UMI COUNT	UMI PERCENT
C:57-G:2579-T:5270	11628	99.61449499	10810	99.60379619
C:57-A:2579-T:5270	16	0.137068449	14	0.128996591
C:57-G:2579-C:5270	15	0.128501671	15	0.138210633
C:57-G:2579-G:5270	6	0.051400668	6	0.055284253
C:57-G:2579-A:5270	3	0.025700334	3	0.027642127

TABLE 15-continued

(4 of 4): P51 variant counts using reads or UMIs at three attenuation sites.				
ATTENUATION SITE	READ COUNT	READ PERCENT	UMI COUNT	UMI PERCENT
T:57-G:2579-T:5270	3	0.025700334	3	0.027642127
A:57-G:2579-T:5270	1	0.008566778	1	0.009214042
C:57-C:2579-T:5270	1	0.008566778	1	0.009214042

[0369] The Pacbio sequencing statistics are shown in Table 16. The P1, P3, P5 and P51 reads were aligned to the dengue virus reference genome (GenBank accession U87411; Dengue Virus type 2, Strain 16681) using BLASR. The numbers of mapped reads are shown in Table 16.

[0370] Similar to the revertant analysis using the BLASR and Quiver workflow, haplotypes T:57-G:2579-T:5270,

T:57-A:2579-A:5270, C:57-A:2579-T:5270 and C:57-G:2579-T:5270 were the most abundant in P1, P3, P5 and P51, respectively. But in contrast to the BLASR and Quiver workflow, which identified 2-5 haplotypes per revertant, the more recent UMI workflow identified 8-11 haplotypes per revertant (see Table 15).

TABLE 16

Pacbio sequencing statistics for mono-revertants P1, P3, P5 and P51.					
Sample	Yield (Mbps)	Reads of Insert	Mean number of Passes	Mean Read Quality	Mean Read Of Insert Length (bps)
P1	81.9	1754	10.15	0.9935	4669
P3	52	11134	10.56	0.9948	4674
P5	66	14072	9.24	0.9924	4699
P51	65.6	12199	8.67	0.9955	5378

[0371] Analysis of the results showed expected changes in nucleotides at three main attenuation loci after RT-PCR and sequencing of these revertants. A total read of 8.7 K post filtration (size 5.4 Kb) using 6500 copies of viral RNA in a 30 cycle PCR reaction has been obtained.

Quantitation of Percentage of Revertants in a Complex Mixture

Preparation of Complex Mixtures:

[0372] qRT-PCR was conducted on individual revertants and 4 separate mixes were prepared based on the genome copies as given in the following chart (Table 17).

TABLE 17

Combinations of complex samples					
#1	P1 35%	P3 50%	P5 5%	TDV-2 10%	Three Single reversions on three different genomes within population mixed with non reverted TDV2 genomes
#2	P3 10%	P51 20%	P513 40%	TDV-2 30%	A mixture of genomes carrying one, two, or three reversions on three different genomes with population mixed with non reverted TDV2 genomes
#3	P3 10%	P51 20%	P513 40%	wt DENV-2 30%	A mixture of genomes carrying one, two, or three reversions on three different genomes within population mixed with wild type DENV-2 genomes (co-infection)
#4	wt DEN-2 75%	TDV-2 25%	A mixture of genomes carrying wild type dengue 2 genome due to secondary (co-infection) infection within population mixed with non reverted TDV2 genomes		
Revertants		Reversion site		Number of Reversion sites	
P1		N51		1	
P3		N53		1	
P513		5'NCR		1	
P51		5'NCR and N51		2	
P513		5'NCR, N51 and N53		3	
TDV-2		No Reversions			
Wt DEN-2		Would be similar to P513 revertant			

[0373] Large Scale PCR to Incorporate UMIs into Complex Revertant Mixtures PCR products were run in a 100 TAE agarose gel. The results show a single amplicon of expected size for Complex Mix 1-4. Large scale RT-PCR to incorporate UMI into complex revertant mixtures, complex mix 1-4. RT was done using the SMARTER IIA Oligonucle-

otide, followed by touchdown PCR with the 5F forward primer and the 5' PCR Primer IIA reverse primer.

Pacbio Sequencing and Variant Analysis

[0374] One SMRTcell was sequenced per complex revertant mixture. The Pacbio sequencing statistics are shown in Table 18.

TABLE 18

Pacbio sequencing stats for four complex revertant pools.					
Sample	Yield (Mbps)	Reads of Insert	Mean number of Passes	Mean Read Quality	Mean Read Of Insert Length (bps)
P#1	52.3	9737	8.57	0.9951	5375
P#2	44.5	8959	9.5	0.9961	4967
P#3	49.7	9254	9.12	0.996	5380
P#4	44.1	8195	8.95	0.9967	5388

Table 19: Variant Counts Tables Using Reads or UMIs at Three Attenuation Loci for Complex Mixtures of TDV2 Mono Revertants

TABLE 19

(1 of 4): Complex Mix 1 (P#1) variant counts using reads or UMIs at three attenuation sites. The estimated % of haplotypes in the mix were 23.27 (P1), 57.76 (P3), 6.89 (P5) and 11.20 (TDV2).				
P#1 ATTENUATION SITE	READ COUNT	READ PERCENT	UMI COUNT	UMI PERCENT
T:57-A:2579-A:5270	5261	57.78143877	4955	57.76404756
T:57-G:2579-T:5270	2117	23.25096101	1996	23.26882723
T:57-A:2579-T:5270	1025	11.2575508	961	11.20307764
C:57-A:2579-T:5270	625	6.864360242	591	6.889717883
T:57-G:2579-A:5270	43	0.472267985	42	0.489624621
C:57-G:2579-T:5270	9	0.098846787	9	0.104919562
C:57-A:2579-A:5270	8	0.087863811	8	0.093261833
T:57-A:2579-G:5270	8	0.087863811	7	0.081604104
C:57-A:2579-C:5270	3	0.032948929	3	0.034973187
G:57-A:2579-A:5270	3	0.032948929	3	0.034973187
A:57-G:2579-T:5270	1	0.010982976	1	0.011657729
T:57-T:2579-A:5270	1	0.010982976	1	0.011657729
T:57-T:2579-T:5270	1	0.010982976	1	0.011657729

TABLE 19

(2 of 4): Complex Mix 2 (P#2) variant counts using reads or UMIs at three attenuation sites. The estimated % of haplotypes in the mix were 10.84 (P3), 21.56 (P51), 37.16 (P513) and 29.59 (TDV-2).				
P#2 ATTENUATION SITE	READ COUNT	READ PERCENT	UMI COUNT	UMI PERCENT
C:57-G:2579-A:5270	2663	37.14086471	2558	37.16402731
T:57-A:2579-T:5270	2119	29.55369596	2037	29.59465349
C:57-G:2579-T:5270	1547	21.57601116	1484	21.56036612
T:57-A:2579-A:5270	780	10.87866109	746	10.83829725
T:57-G:2579-T:5270	17	0.237099024	16	0.232456778
T:57-G:2579-A:5270	11	0.153417015	11	0.159814035
C:57-A:2579-T:5270	9	0.125523013	9	0.130756937
C:57-A:2579-A:5270	8	0.111576011	6	0.087171292
T:57-A:2579-G:5270	5	0.069735007	5	0.072642743
C:57-G:2579-C:5270	2	0.027894003	2	0.029057097
C:57-T:2579-T:5270	2	0.027894003	2	0.029057097
A:57-A:2579-T:5270	1	0.013947001	1	0.014528549
C:57-G:2579-G:5270	1	0.013947001	1	0.014528549
G:57-A:2579-A:5270	1	0.013947001	1	0.014528549
G:57-A:2579-T:5270	1	0.013947001	1	0.014528549

TABLE 19-continued

(2 of 4): Complex Mix 2 (P#2) variant counts using reads or UMIs at three attenuation sites. The estimated % of haplotypes in the mix were 10.84 (P3), 21.56 (P51), 37.16 (P513) and 29.59 (TDV-2).				
P#2	READ			UMI
ATTENUATION SITE	COUNT	READ PERCENT	UMI COUNT	PERCENT
T:57-A:2579-C:5270	1	0.013947001	1	0.014528549
T:57-C:2579-T:5270	1	0.013947001	1	0.014528549
T:57-T:2579-T:5270	1	0.013947001	1	0.014528549

TABLE 19

(3 of 4): Complex Mix 3 (P#3) variant counts using reads or UMIs at three attenuation sites. The estimated % of haplotypes in the mix were 9.62 (P3), 25.12 (P51), 37.16 (P513) and 27.72 (wt DENV-2).				
P#3	READ			UMI
ATTENUATION SITE	COUNT	READ PERCENT	UMI COUNT	PERCENT
C:57-G:2579-A:5270	5740	64.75631769	5428	64.8894202
C:57-G:2579-T:5270	2237	25.23691336	2101	25.11655708
T:57-A:2579-A:5270	853	9.623194946	805	9.623430962
C:57-A:2579-A:5270	13	0.14666065	12	0.143454871
C:57-G:2579-C:5270	5	0.056407942	5	0.059772863
C:57-G:2579-G:5270	5	0.056407942	4	0.04781829
T:57-G:2579-A:5270	4	0.045126354	4	0.04781829
A:57-G:2579-A:5270	2	0.022563177	2	0.023909145
T:57-A:2579-T:5270	2	0.022563177	2	0.023909145
T:57-G:2579-T:5270	2	0.022563177	1	0.011954573
C:57-A:2579-T:5270	1	0.011281588	1	0.011954573

TABLE 19

(4 of 4): Complex Mix 4 (P#4) variant counts using reads or UMIs at three attenuation sites. The estimated % of haplotypes in the mix were 65.57 (wt DENV-2) and 33.70 (TDV-2).				
P#4	READ			UMI
ATTENUATION SITE	COUNT	READ PERCENT	UMI COUNT	PERCENT
C:57-G:2579-A:5270	5159	65.64448403	4926	65.57507987
T:57-A:2579-T:5270	2645	33.65568138	2532	33.70607029
C:57-A:2579-T:5270	11	0.139966917	11	0.146432375
C:57-G:2579-T:5270	9	0.114518387	9	0.119808307
T:57-G:2579-A:5270	8	0.101794121	8	0.106496273
C:57-A:2579-A:5270	6	0.076345591	6	0.079872204
T:57-A:2579-A:5270	5	0.063621326	5	0.06656017
C:57-G:2579-G:5270	4	0.050897061	4	0.053248136
T:57-A:2579-C:5270	4	0.050897061	3	0.039936102
T:57-G:2579-T:5270	3	0.038172796	3	0.039936102
T:57-A:2579-G:5270	2	0.02544853	2	0.026624068
A:57-G:2579-A:5270	1	0.012724265	1	0.013312034
C:57-T:2579-A:5270	1	0.012724265	1	0.013312034
T:57-C:2579-T:5270	1	0.012724265	1	0.013312034

[0375] For Complex Mix 1, 13 haplotypes were identified with T:57-A:2579-A:5270 being the most abundant, followed by T:57-G:2579-T:5270, T:57-A:2579-T:5270 and C:57-A:2579-T:5270 (Table 19; 1 of 4); the remaining haplotypes made up a very small percentage of the population, ~0.77%. Complex Mix 2 contained the largest number of haplotypes, 18, with C:57-G:2579-A:5270 being the most abundant followed by T:57-A:2579-T:5270, C:57-G:2579-T:5270 and T:57-A:2579-A:5270; the remaining haplotypes made up ~0.85% of the population (Table 19; 2 of 4).

Complex Mix 3 contained the fewest numbers of haplotypes with 11. Haplotype C:57-G:2579-A:5270 was the most abundant, followed by C:57-G:2579-T:5270, and T:57-A:2579-A:5270, with the remaining haplotype making up ~0.38% of the population (Table 19; 3 of 4). Complex Mix 4 contained the second most haplotypes with 14. However, just two haplotypes, C:57-G:2579-A:5270 and T:57-A:2579-T:5270, made up over 99% of the population with the remaining haplotypes making up only ~0.7% of the population (Table 19; 4 of 4).

TABLE 20

Read and UMI counting at attenuation sites for complex samples				
ATTENUATION SITES	READ_COUNT	READ_PERCENT	UMI_COUNT	UMI_PERCENT
P#1				
T:57-A:2579-A:5270	5261	57.78	4955	57.76
T:57-G:2579-T:5270	2117	23.25	1996	23.26
T:57-A:2579-T:5270	1025	11.25	961	11.20
C:57-A:2579-T:5270	625	6.86	891	6.88
P#2				
C:57-G:2579-A:5270	2663	37.14	2558	37.16
T:57-A:2579-T:5270	2119	29.55	2037	29.59
C:57-G:2579-T:5270	1847	21.57	1484	21.56
T:57-A:2579-A:5270	780	10.87	746	10.83
P#3				
C:57-G:2579-A:5270	5740	64.75	5428	84.88
C:57-G:2579-T:5270	2237	26.23	2101	25.11
T:57-A:2579-A:5270	853	9.62	805	9.62
P#4				
C:57-G:2579-A:5270	5159	65.64	4926	65.57
T:57-A:2579-T:5270	2645	33.66	2532	33.70

Wild-Type: C:57-G:2579-A:5270

Vaccine: T:57-A:2579-T:5270

Comparison of Pacbio Results with Blinded Complex Samples

[0376] The percentage of each revertant in the complex mix determined by the pacbio sequencing was compared with the blinded mix prepared by Takeda. As can be seen from the Table 21, the quantity of each haplotype in the pools P #1, P #2, P #3 and P #4 was comparable to that of original amounts formulated by Takeda.

TABLE 21

Quantitation of haplotypes in the complex pools by pacbio sequencing				
Pool #	Components of mix			
1	P1	P3	P5	TDV-2
	35	50	5	10
	23.27	57.76	6.89	11.20
2	P3	P51	P513	TDV-2
	10	20	40	30
	10.84	21.58	37.16	29.59
3	P3	P51	P513	Wt DENV-2
	10	20	40	30
	9.62	25.12	64.89	(37.16 + 27.72)
4	Wt DENV-2	TDV-2	—	—
	75	25	—	—
	65.57	33.70	—	—

Identity by RT-PCR Amplification of E Gene Sequences

[0377] Reverse transcription polymerase chain reaction (RT-PCR) is a sensitive method for the detection of RNA, which involves two step reactions. RNA is first reverse transcribed to cDNA by reverse transcriptase. The cDNA is then used as template for the subsequent PCR step. Amplified DNA fragments are separated and visualized in a 2% agarose gel; method conditions are listed in Table 23. The identity of each monovalent drug substance is determined by a positive result for RT-PCR amplification of an E gene sequence specific to the appropriate dengue serotype, and a negative result for RT-PCR amplification of the E gene sequences specific to the other three dengue serotypes. Four primer sets were designed to differentiate the four serotypes of the dengue vaccine viruses. The serotype-specific primer sets were designed such that each primer has only one specific binding site within the genome of the corresponding dengue serotype and each primer set generates a DNA fragment of specific size. Viral RNA is isolated from the test sample using the QIAamp® cador® Pathogen Kit according to the manufacturer's instructions. RT-PCR on the isolated viral RNA is performed using a one-step RT-PCR kit (highQu GmbH). The amplicons for the four vaccine virus strains are 419, 390, 731 and 226 base pairs for TDV-1, TDV-2, TDV-3 and TDV-4, respectively, and are easily differentiated on a 2% agarose gel (Table 22).

TABLE 22

Primer Sequences and Size of Amplicons for RT-PCR-based Identity Test			
Virus Strain		Serotype-Specific E Gene Primer Sequences	Amplicon Size (in Basepairs)
TDV-1	Forward	5'-CCGACTACGGAACCTTACA-3'	419
	Reverse	5'-GGTCTCAGCCACTTCTTCTCTA-3'	
TDV-2	Forward	5'-GCCATGCACACAGCACTTAC-3'	390
	Reverse	5'-CTTCAGTTGTCCCGGCTCTA-3'	

TABLE 22-continued

Primer Sequences and Size of Amplicons for RT-PCR-based Identity Test			
Virus Strain		Serotype-Specific E Gene Primer Sequences	Amplicon Size (in Basepairs)
TDV-3	Forward	5'-AAGCCACGCTGGACATAG-3'	731
	Reverse	5'-AAGTCCCCGCAAGATACT-3'	
TDV-4	Forward	5'-GCCGGAGAGATGTGGTAGAC-3'	226
	Reverse	5'-CTGGGAGTTATCGTGGCTGTA-3'	

[0378] The above primer sequences have been assigned the following sequence identifiers.

TDV-1	Forward	SEQ ID NO: 24
	Reverse	SEQ ID NO: 25
TDV-2	Forward	SEQ ID NO: 26
	Reverse	SEQ ID NO: 27
TDV-3	Forward	SEQ ID NO: 28
	Reverse	SEQ ID NO: 29
TDV-4	Forward	SEQ ID NO: 30
	Reverse	SEQ ID NO: 31

TABLE 23

RT-PCR Condition		
Temperature	Time	Cycles
50° C.	20 min	1
95° C.	2 min	1
95° C.	30 sec	35
66° C.	30 sec	
72° C.	1 min	
4° C.	0 to infinity	1

Residual Host Cell DNA by qPCR

[0379] A qPCR based method is used for residual host cell DNA that detects and estimates Vero DNA. This method utilizes oligo-magnetic bead capture to improve DNA detection and is available in a kit manufactured by Life Technologies Inc. The method was employed to ensure low levels of residual Vero host cell DNA can be detected in the bulk drug substance.

[0380] For DNA isolation and purification, the PrepSEQ® Residual DNA Sample Preparation Kit extracts host cell DNA from products produced in Vero cells. The kit uses chemical lysis and magnetic beads to efficiently extract genomic DNA from diverse sample types, including samples

that contain high protein and low DNA concentration. Each test sample is extracted in triplicate and a single PCR reaction is performed for each extraction.

[0381] For quantitative PCR, the resDNASEQ® Quantitative DNA kit is used to quantify residual DNA from the Vero cell line which is used for production of bulk drug substance. The resDNASEQ® Quantitative DNA kit uses TaqMan® quantitative PCR to perform rapid, specific quantitation of low levels of residual host cell DNA in testing samples by interpolation of the standard curve provided in the kit.

[0382] There are 4 controls used for the DNA extraction and qPCR. Extraction Recovery Control (ERC) is drug substance sample spiked with 500 µg of Vero Control DNA. Negative Control (NC) is molecular biology grade water used as negative control for the whole procedure, from DNA extraction to qPCR. No Template Control (NTC) is molecular biology grade water as negative control in qPCR. Internal Positive Control (IPC) is provided in the qPCR kit as positive control for qPCR.

[0383] For the testing to be valid, the standard curve generated must fulfill the conditions of R² more than 0.99, efficiency of 90% to 110%, and y-intercept between 24 and 29. The NTC and SE must have a mean of triplicate negative controls DNA quantity of less than 30 fg. The ERC must have percentage recovery mean of triplicate between 50% and 150%. The IPC must have a Ct value less than 40.00 in all test samples and controls. For the sample result to be valid, the repeatability of the test sample replicates must be less than or equal to 25% CV.

[0384] The Ct values of valid test samples from the PCR are used to determine the result. The result is reported as 'Negative' if the Ct is undetermined. The result is reported as 'Positive' if the Ct values of at least 2 out of 3 replicates fall within the standard curve. The DNA concentration will be calculated by interpolating the Ct values from the standard curve.

SEQUENCE LISTING

```

<160> NUMBER OF SEQ ID NOS: 31

<210> SEQ ID NO 1
<211> LENGTH: 24
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 1

aagacagatt ctttgaggga gcta

```

-continued

<210> SEQ ID NO 2
<211> LENGTH: 20
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 2

ggtctgtgaa atgggcttcg 20

<210> SEQ ID NO 3
<211> LENGTH: 23
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 3

gacagattct ttgagggagc taa 23

<210> SEQ ID NO 4
<211> LENGTH: 20
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 4

gggtctgtga aatgggcttc 20

<210> SEQ ID NO 5
<211> LENGTH: 23
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 5

acagattctt tgagggagct aag 23

<210> SEQ ID NO 6
<211> LENGTH: 21
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 6

gtctgtgaaa tgggcttcgt c 21

<210> SEQ ID NO 7
<211> LENGTH: 25
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 7

agacagattc tttgagggag ctaag 25

<210> SEQ ID NO 8
<211> LENGTH: 20

-continued

<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 8

gtctgtgaaa tgggcttcgt 20

<210> SEQ ID NO 9
<211> LENGTH: 24
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 9

gacagattct ttgaggagc taag 24

<210> SEQ ID NO 10
<211> LENGTH: 21
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 10

cttgctgggt ctgtgaaatg g 21

<210> SEQ ID NO 11
<211> LENGTH: 25
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 11

aagacagatt ctttgaggga gctaa 25

<210> SEQ ID NO 12
<211> LENGTH: 21
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 12

tgggtctgtg aaatgggctt c 21

<210> SEQ ID NO 13
<211> LENGTH: 24
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 13

agttgttagt ctacgtggac cgac 24

<210> SEQ ID NO 14
<211> LENGTH: 30
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

-continued

<400> SEQUENCE: 14

gtgacactgt atttgggagt catggtgcag

30

<210> SEQ ID NO 15

<211> LENGTH: 25

<212> TYPE: DNA

<213> ORGANISM: artificial

<220> FEATURE:

<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 15

ggtaatggtg ttgttacaag gagtg

25

<210> SEQ ID NO 16

<211> LENGTH: 24

<212> TYPE: DNA

<213> ORGANISM: artificial

<220> FEATURE:

<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 16

aagcccttag aggacttcca ataa

24

<210> SEQ ID NO 17

<211> LENGTH: 24

<212> TYPE: DNA

<213> ORGANISM: artificial

<220> FEATURE:

<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 17

tcagcatcct tccaatctct ttcc

24

<210> SEQ ID NO 18

<211> LENGTH: 24

<212> TYPE: DNA

<213> ORGANISM: artificial

<220> FEATURE:

<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 18

ccgcagagat cgttttcctg cctg

24

<210> SEQ ID NO 19

<211> LENGTH: 25

<212> TYPE: DNA

<213> ORGANISM: artificial

<220> FEATURE:

<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 19

ctaactggtg atagcagcct catgg

25

<210> SEQ ID NO 20

<211> LENGTH: 28

<212> TYPE: DNA

<213> ORGANISM: artificial

<220> FEATURE:

<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 20

gaaatgggct tcgtccatga taatcagg

28

-continued

<210> SEQ ID NO 21
 <211> LENGTH: 25
 <212> TYPE: DNA
 <213> ORGANISM: artificial
 <220> FEATURE:
 <223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 21

aagcagtggg atcaacgcag agtac 25

<210> SEQ ID NO 22
 <211> LENGTH: 62
 <212> TYPE: DNA
 <213> ORGANISM: artificial
 <220> FEATURE:
 <223> OTHER INFORMATION: synthetic
 <220> FEATURE:
 <221> NAME/KEY: misc_feature
 <222> LOCATION: (26)..(41)
 <223> OTHER INFORMATION: n is a, c, g, or t

<400> SEQUENCE: 22

aagcagtggg atcaacgcag agtacnnnnn nnnnnnnnnn ngctctgtgaa atgggcttcg 60

tc 62

<210> SEQ ID NO 23
 <211> LENGTH: 10723
 <212> TYPE: DNA
 <213> ORGANISM: dengue virus type 2

<400> SEQUENCE: 23

agttgttagt ctactgggac cgacaaagac agattctttg agggagctaa gctcaatgta 60

gttctaacag ttttttaatt agagagcaga tctctgatga ataaccaacg gaaaaaggcg 120

aaaaaacgcg ctttcaatat gctgaaacgc gagagaaacc gcgtgtcgac tgtgcaacag 180

ctgacaaaga gattctcact tggaatgctg cagggacgag gaccattaaa actgttcacg 240

gccctgggtg cgcttcctcg tttcctaaca atcccaccaa cagcagggat attgaagaga 300

tggggaacaa ttaaaaaatc aaaagctatt aatgttttga gaggggttcag gaaagagatt 360

ggaaggatgc tgaacatctt gaataggaga cgcagatctg caggcatgat cattatgctg 420

attccaacag tgatggcggt ccatttaacc acacgtaacg gagaaccaca catgatcgtc 480

agcagacaag agaaagggaa aagtcttctg tttaaaacag aggttggcgt gaacatgtgt 540

accctcatgg ccatggacct tggatgaattg tgtgaagaca caatcacgta cgagtgtccc 600

cttctcagcg agaatgagcc agaagacata gactgttggt gcaactctac gtccacgtgg 660

gtaacttatg ggacgtgtac caccatggga gaacatagaa gagaaaaaag atcagtggca 720

ctcgttccac atgtgggaat gggactggag acacgaactg aaacatggat gtcacagaa 780

ggggcctgga aacatgtcca gagaattgaa acttggtatc tgagacatcc aggttcacc 840

atgatggcag caatcctggc atacaccata ggaacgacac atttccaaag agcctgatc 900

ttcatcttac tgacagctgt cactccttca atgacaatgc gttgcatagg aatgtcaaat 960

agagactttg tggaagggtt ttcaggagga agctgggttg acatagtctt agaacatgga 1020

agctgtgtga cgacgatggc aaaaaacaaa ccaacattgg attttgaact gataaaaaa 1080

gaagccaaac agcctgccac cctaagggaag tactgtatag aggcaaagct aaccaacaca 1140

-continued

acaacagaat	ctcgctgccc	aacacaaggg	gaacccagcc	taaatgaaga	gcaggacaaa	1200
aggttcgtct	gcaaacactc	catggtagac	agaggatggg	gaaatggatg	tggactattt	1260
ggaaagggag	gcatttgtac	ctgtgtctatg	ttcagatgca	aaaagaacat	ggaaggaaaa	1320
gttgtgcaac	cagaaaaactt	ggaatacacc	atttgtataa	cacctcactc	aggggaagag	1380
catgcagtcg	gaaatgacac	aggaaaacat	ggcaaggaaa	tcaaaataac	accacagagt	1440
tccatcacag	aagcagaatt	gacaggttat	ggcactgtca	caatggagtg	ctctccaaga	1500
acgggcctcg	acttcaatga	gatgggtgtg	ctgcagatgg	aaaataaagc	ttggctggtg	1560
cacaggcaat	ggttcctaga	cctgccgtta	ccatggttgc	ccggagcggg	cacacaaggg	1620
tcaaattgga	tacagaaaga	gacattggtc	actttcaaaa	atccccatgc	gaagaaacag	1680
gatgttggtg	ttttaggatc	ccaagaaggg	gccatgcaca	cagcacttac	aggggccaca	1740
gaaatccaaa	tgtcatcagg	aaacttactc	ttcacaggac	atctcaagtg	caggctgaga	1800
atggacaagc	tacagctcaa	aggaatgtca	tactctatgt	gcacaggaaa	gtttaaagt	1860
gtgaaggaaa	tagcagaaac	acaacatgga	acaatagtta	tcagagtgcg	atatgaaggg	1920
gacggctctc	catgcaagat	cccttttgag	ataatggatt	tggaaaaaag	acatgtctta	1980
ggtcgccctg	ttacagtcaa	cccaattgtg	acagaaaaag	atagcccagt	caacatagaa	2040
gcagaacctc	catttgagga	cagctacatc	atcataggag	tagagccggg	acaactgaag	2100
ctcaactggt	ttaagaaagg	aagttctatc	ggccaaatgt	ttgagacaac	aatgaggggg	2160
gcgaagagaa	tggccatttt	aggtgacaca	gcctgggatt	ttggatcctt	gggaggagtg	2220
tttcatctta	taggaaaggc	tctccaccaa	gtctttggag	caatctatgg	agctgccttc	2280
agtgggggtt	catggactat	gaaaatcctc	ataggagtca	ttatcacatg	gataggaatg	2340
aattcacgca	gcacctcact	gtctgtgaca	ctagtattgg	tgggaattgt	gacactgtat	2400
ttgggagtca	tggtgcaggc	cgatagtggg	tgcgttgtga	gctggaaaaa	caaagaactg	2460
aaatgtggca	gtgggatttt	catcacagac	aacgtgcaca	catggacaga	acaatacaag	2520
ttccaaccag	aatccccctc	aaaactagct	tcagctatcc	agaaagccca	tgaagaggac	2580
atltgtggaa	tccgctcagt	aacaagactg	gagaatctga	tgtggaaaca	aataacacca	2640
gaattgaatc	acattctatc	agaaaatgag	gtgaagttaa	ctattatgac	aggagacatc	2700
aaaggaatca	tgcaggcagg	aaaacgatct	ctgcggcctc	agcccaactg	gctgaagtat	2760
tcattggaaa	catggggcaa	agcaaaaatg	ctctctacag	agtctcataa	ccagaccttt	2820
ctcattgatg	gccccgaaac	agcagaatgc	cccaacacaa	atagagcttg	gaattcgttg	2880
gaagtgaag	actatggctt	tggagtattc	accaccaata	tatggctaaa	attgaaagaa	2940
aaacaggatg	tattctgcga	ctcaaaactc	atgtcagcgg	ccataaaaga	caacagagcc	3000
gtccatgccg	atatgggtta	ttggatagaa	agtgactca	atgacacatg	gaagatagag	3060
aaagcctctt	tcattgaagt	taaaaactgc	cactggccaa	aatcacacac	cctctggagc	3120
aatggagtgc	tagaaagtga	gatgataatt	ccaaagaatc	tcgctggacc	agtgtctcaa	3180
cacaactata	gaccaggcta	ccatacacaa	ataacaggac	catggcatct	aggtaagctt	3240
gagatggact	ttgatttctg	tgatggaaca	acagtggtag	tgactgagga	ctgcggaaat	3300
agaggacctt	ctttgagaac	aacctatgcc	tctggaaaac	tcataacaga	atgggtgctgc	3360
cgatcttgca	cattaccacc	gctaagatac	agaggtgagg	atgggtgctg	gtacgggatg	3420

-continued

gaaatcagac cattgaagga gaaagaagag aatttgggtca actccttggt cacagctgga	3480
catgggcagg tegacaactt ttcactagga gtcttgggaa tggcattggt cctggaggaa	3540
atgcttagga cccgagtagg aacgaaacat gcaatactac tagttgcagt ttcttttgtg	3600
acattgatca cagggaacat gtcccttaga gacctgggaa gagtgatggt tatggtaggc	3660
gccactatga cggatgacat aggtatgggc gtgacttata ttgccctact agcagccttc	3720
aaagtcagac caacttttgc agctggacta ctcttgagaa agctgacctc caaggaattg	3780
atgatgacta ctataggaat tgtactctc tcccagagca ccataccaga gaccattctt	3840
gagttgactg atgcgttagc cttaggcatg atggctctca aaatggtgag aaatatggaa	3900
aagtatcaat tggcagtgac tatcatggct atcttggtgcg tcccaaacgc agtgatatta	3960
caaaacgcac ggaaagttag ttgcacaata ttggcagtggt tgtccgttcc cccactgttc	4020
ttaacatcct cacagcaaaa aacagattgg ataccattag cattgacgat caaaggtctc	4080
aatccaacag ctatttttct aacaaccctc tcaagaacca gcaagaaaag gagctggcca	4140
ttaaagtagg ctatcatggc agtcgggatg gtgagcattt tagccagtcc tctcctaaaa	4200
aatgatattc ccatgacagg accatttagt gctggagggc tctcactgt gtgctacgtg	4260
ctcactggac gatcgccga tttggaactg gagagagcag ccgatgtcaa atgggaagac	4320
caggcagaga tatcaggaag cagtccaatc ctgtcaataa caatatcaga agatggtagc	4380
atgtcgataa aaaatgaaga ggaagaacaa aactgacca tactcattag aacaggattg	4440
ctggtgatct caggactttt tctgtatca ataccaatca cggcagcagc atggtacctg	4500
tgggaagtga agaaacaacg ggcgggagta ttgtgggatg ttcttcacc cccacccatg	4560
ggaaaggctg aactggaaga tggagcctat agaattaagc aaaaagggat tcttgatat	4620
tcccagatcg gagccggagt ttacaaagaa ggaacattcc atacaatgtg gcatgtcaca	4680
cgtggcgctg ttctaataca taaaggaaaag aggattgaac catcatgggc ggacgtcaag	4740
aaagacctaa tatcatatgg aggagctgg aagttagaag gagaatggaa ggaaggagaa	4800
gaagtccagg tattggcact ggagcctgga aaaaatccaa gagccgtcca aacgaaacct	4860
ggtcttttca aaaccaacgc cggacaataa ggtgctgtat ctctggactt ttctcctgga	4920
acgtcaggat ctccaattat cgacaaaaaa ggaaaagttg tgggtcttta tggtaatggt	4980
gttgttacaa ggagtggagc atatgtgagt gctatagccc agactgaaaa aagcattgaa	5040
gacaaccag agatcgaaga tgacatttcc cgaaagagaa gactgaccat catggacctc	5100
caccagagag cgggaaagac gaagagatac ctcccgcca tagtcagaga agctataaaa	5160
cggggtttga gaacattaat cttggcccc actagagttg tggcagctga aatggaggaa	5220
gcccttagag gacttccaat aagataccag acccagcca tcagagctgt gcacaccggg	5280
cgggagattg tggacctaat gtgtcatgcc acatttacca tgaggctgct atcaccagtt	5340
agagtgccaa actacaacct gattatcatg gacgaagccc atttcacaga cccagcaagt	5400
atagcagcta gaggatacat ctcaactcga gtggagatgg gtgaggcagc tgggattttt	5460
atgacagcca cccccggg aagcagagac ccatttcctc agagcaatgc accaatcata	5520
gatgaagaaa gagaaatccc tgaacgctcg tggaaatccg gacatgaatg ggtcacggat	5580
tttaaaggga agactgtttg gttcgttcca agtataaaa caggaaatga tatagcagct	5640
tgctgagga aaaatggaaa gaaagtgata caactcagta ggaagacctt tgattctgag	5700

-continued

tatgtcaaga	ctagaaccaa	tgattgggac	ttcgtgggta	caactgacat	ttcagaaatg	5760
gggtccaatt	tcaaggctga	gaggggtata	gaccccagac	gctgcatgaa	accagtcata	5820
ctaacagatg	gtgaagagcg	gggtattctg	gcaggaccta	tgccagtgc	ccactctagt	5880
gcagcacaaa	gaagaggag	aataggaaga	aatccaaaaa	atgagaatga	ccagtacata	5940
tacatggggg	aacctctgga	aatgatgaa	gactgtgcac	actggaaaga	agctaaaatg	6000
ctcctagata	acatcaacac	gccagaagga	atcattccta	gcatgttcga	accagagcgt	6060
gaaaagggtg	atgccattga	tggcgaatac	cgcttgagag	gagaagcaag	gaaaaccttt	6120
gtagacttaa	tgagaagagg	agacctacca	gtctgggttg	cctacagagt	ggcagctgaa	6180
ggcatcaact	acgcagacag	aagggtggtg	tttgatggag	tcaagaacaa	ccaaatccta	6240
gaagaaaacg	tggaagtga	aatctggaca	aaagaagggg	aaaggaagaa	attgaaaccc	6300
agatgggttg	atgctaggat	ctattctgac	ccactggcgc	taaaagaatt	taaggaattt	6360
gcagccggaa	gaaagtctct	gacctgaac	ctaatacag	aaatgggtag	gctcccaacc	6420
ttcatgactc	agaaggcaag	agacgcactg	gacaacttag	cagtgtctga	cacggctgag	6480
gcaggtggaa	gggcgtacaa	ccatgctctc	agtgaactgc	cggagaccct	ggagacattg	6540
cttttactga	cacttctggc	tacagtcacg	ggagggatct	ttttattctt	gatgagcgca	6600
aggggcatag	ggaagatgac	cctgggaatg	tgctgcataa	tcacggctag	catcctccta	6660
tggtacgcac	aaatacagcc	acactggata	gcagcttcaa	taatactgga	gttttttctc	6720
atagttttgc	ttattccaga	acctgaaaaa	cagagaacac	cccaagacaa	ccaaactgacc	6780
tacgttgtca	tagccatcct	cacagtgggtg	gccgcaacca	tggcaaacga	gatgggtttc	6840
ctagaaaaaa	cgaagaaaga	tctcggattg	ggaagcattg	caaccagca	accagagagc	6900
aacatcctgg	acatagatct	acgtcctgca	tcagcatgga	cgctgtatgc	cgtggccaca	6960
acatttgtta	caccaatggt	gagacatagc	attgaaaatt	cctcagtga	tgtgtcccta	7020
acagctatag	ccaaccaagc	cacagtgtta	atgggtctcg	ggaaaggatg	gccattgtca	7080
aagatggaca	tccgagttcc	ccttctcgcc	attggatgct	actcacaagt	caaccccata	7140
actctcacag	cagctctttt	cttattggta	gcacattatg	ccatcatagg	gccaggactc	7200
caagcaaaag	caaccagaga	agctcagaaa	agagcagcgg	cgggcacatc	gaaaaacca	7260
actgtcgatg	gaataaacgt	gattgacctc	gatccaatac	cttatgatcc	aaagtttgaa	7320
aagcagttgg	gacaagtaat	gctcctagtc	ctctgcgtga	ctcaagtatt	gatgatgagg	7380
actacatggg	ctctgtgtga	ggctttaacc	ttagctaccg	ggcccatctc	cacattgtgg	7440
gaaggaaatc	cagggaggtt	ttggaacact	accattgctg	tgtcaatggc	taacattttt	7500
agagggagtt	acttgcccg	agctggactt	ctcttttcta	ttatgaagaa	cacaaccaac	7560
acaagaaggg	gaactggcaa	cataggagag	acgcttgag	agaaatggaa	aagccgattg	7620
aacgcattgg	gaaaaagtga	attccagatc	tacaagaaaa	gtggaatcca	ggaagtggat	7680
agaaccttag	caaaagaagg	cattaaaaga	ggagaaacgg	accatcacgc	tgtgtcgcg	7740
ggctcagcaa	aactgagatg	gttcgttgag	agaaacatgg	tcacaccaga	agggaaagta	7800
gtggacctcg	gttggtggcag	aggaggctgg	tcatactatt	gtggaggact	aaagaatgta	7860
agagaagtca	aaggcctaac	aaaaggagga	ccaggacacg	aagaacccat	ccccatgtca	7920
acatatgggt	ggaatctagt	gcgtcttcaa	agtgaggttg	acgttttctt	catcccgcca	7980

-continued

gaaaagtgtg acacattatt gtgtgacata ggggagtcac caccaaatcc cacagtggaa	8040
gcaggacgaa cactcagagt ccttaactta gtagaaaatt ggttgaacaa caaactcaa	8100
ttttgcataa aggttctcaa cccatatatg ccctcagtc tagaaaaat ggaagcacta	8160
caaaggaaat atggaggagc cttagttagg aatccactct cagaaaactc cacacatgag	8220
atgtactggg tatccaatgc ttccgggaac atagtgtcat cagtgaacat gatttcaagg	8280
atgttgatca acagatttac aatgagatac aagaagcca cttacgagcc ggatgttgac	8340
ctcggaagcg gaaccgtaa catcgggatt gaaagtgaga taccaaacct agatataatt	8400
gggaaaagaa tagaaaaat aaagcaagag catgaaacat catggcacta tgaccaagac	8460
caccataca aaacgtgggc ataccatggc agctatgaaa caaacagac tggatcagca	8520
tcatccatgg tcaacggagt ggtcaggctg ctgacaaaac cttgggacgt cgtcccatg	8580
gtgacacaga tggcaatgac agacacgact ccatttggac aacagcgcgt ttttaagag	8640
aaagtggaca cgagaacca agaaccgaaa gaaggcacga agaaactaat gaaaataaca	8700
gcagagtggc tttggaaaga attagggag aaaaagacac ccaggatgtg caccagagaa	8760
gaattcaca gaaaggtgag aagcaatgca gccttggggg ccgtattcac tgatgagaa	8820
aagtggaaat cggcacgtga ggtgttgaa gatagtaggt tttgggagct ggttgacaag	8880
gaaaggaatc tccatcttga aggaaagtgt gaaacatgtg tgtacaacat gatgggaaaa	8940
agagagaaga agctagggga attcggcaag gcaaaaggca gcagagccat atggtacatg	9000
tggcttggag cacgcttctt agagtttgaa gccctaggat tcttaaatga agatcactgg	9060
ttctccagag agaactccct gagtggagtg gaaggagaag ggctgcacaa gctaggttac	9120
attctaagag acgtgagcaa gaaagaggga ggagcaatgt atgccgatga caccgcagga	9180
tgggatacaa gaatcacact agaagaccta aaaaatgaag aaatggtaac aaaccacatg	9240
gaaggagaac acaagaaact agccgaggcc attttcaaac taacgtacca aaacaagggtg	9300
gtgcgtgtgc aaagaccaac accaagaggc acagtaatgg acatcatatc gagaagagac	9360
caaagaggta gtggacaagt tggcaccat ggactcaata ctttcaccaa tatggaagcc	9420
caactaatca gacagatgga gggagaagga gtctttaaaa gcattcagca cctaacaatc	9480
acagaagaaa tcgctgtgca aaactgggta gcaagagtgg ggcgcgaaag gttatcaaga	9540
atggccatca gtggagatga ttgtgtgtg aaacctttag atgacagggt cgcaagcgct	9600
ttaacagctc taaatgacat gggaaagatt aggaaagaca tacaacaatg ggaaccttca	9660
agaggatgga atgattggac acaagtgcc ttctgttcac accatttcca tgagttaatc	9720
atgaaagacg gtcgcgtact cgttgttcca tgtagaaacc aagatgaact gattggcaga	9780
gcccgaatct cccaaggagc aggggtggct ttgcgggaga cggcctgttt ggggaagtct	9840
tacgccccaa tgtggagctt gatgtacttc cacagacgcg acctcaggct ggcggcaaat	9900
gctatttgc tggcagttac atcacattgg gttccaacaa gtcgaacaac ctggtccata	9960
catgctaacc atgaatggat gacaacggaa gacatgctga cagtctgga cagggtgtgg	10020
attcaagaaa acccatggat ggaagacaaa actccagtgg aatcatggga ggaaatccca	10080
tacttgggga aaagagaaga ccaatgggtc ggctcattga ttgggttaac aagcagggcc	10140
acctgggcaa agaacatcca agcagcaata aatcaagtta gatcccttat aggcattgaa	10200
gaatacacag attacatgcc atccatgaaa agattcagaa gagaagagga agaagcagga	10260

-continued

gtttctgtggt agaaagcaaa actaacatga aacaaggcta gaagtcaggt cggattaagc 10320
catagtacgg aaaaaactat gctacctgtg agcccccgtcc aaggacgtta aaagaagtca 10380
ggccatcata aatgccatag cttgagtaaa ctatgcagcc ttagctccca cctgagaagg 10440
tgtaaaaaat cggggaggcc acaaaccatg gaagctgtac gcatggcgta gtggactagc 10500
ggtttagagga gaccctctcc ttacaaatcg cagcaacaat gggggcccaa ggcgagatga 10560
agctgtagtgc tcgctggaag gactagaggt tagaggagac ccccccgaac caaaaaacag 10620
catattgacg ctgggaaaga ccagagatcc tgctgtctcc tcagcatcat tccaggcaca 10680
gaacgccaga aaatggaatg gtgctgttga atcaacaggt tct 10723

<210> SEQ ID NO 24
<211> LENGTH: 20
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 24

ccgactacgg aacccttaca 20

<210> SEQ ID NO 25
<211> LENGTH: 23
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 25

ggtctcagcc acttctttct cta 23

<210> SEQ ID NO 26
<211> LENGTH: 22
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 26

gccatgcaca cacagcactt ac 22

<210> SEQ ID NO 27
<211> LENGTH: 20
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 27

cttcagttgt cccggtctca 20

<210> SEQ ID NO 28
<211> LENGTH: 19
<212> TYPE: DNA
<213> ORGANISM: artificial
<220> FEATURE:
<223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 28

aagcccacgc tggacatag 19

-continued

<210> SEQ ID NO 29
 <211> LENGTH: 20
 <212> TYPE: DNA
 <213> ORGANISM: artificial
 <220> FEATURE:
 <223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 29

aagtgcccg caaagatact

20

<210> SEQ ID NO 30
 <211> LENGTH: 20
 <212> TYPE: DNA
 <213> ORGANISM: artificial
 <220> FEATURE:
 <223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 30

gccggagaga tgtggtagac

20

<210> SEQ ID NO 31
 <211> LENGTH: 21
 <212> TYPE: DNA
 <213> ORGANISM: artificial
 <220> FEATURE:
 <223> OTHER INFORMATION: synthetic

<400> SEQUENCE: 31

ctgggagtta tcgtggctgta

21

1. A method for detecting the presence of at least one revertant virus in a sample containing a live, attenuated virus comprising the steps:

- a) purifying virus nucleic acids from said sample;
- b) optionally, preparing double-stranded DNA templates from said virus nucleic acids;
- c) amplifying said viral nucleic acids of a) or said double-stranded DNA templates of b) by polymerase chain reaction (PCR) thereby generating PCR products having a length of at least 2 kb and comprising at least one attenuation locus of said attenuated virus;
- d) sequencing said PCR products by single-molecule real-time (SMRT) sequencing; and
- e) comparing at least one nucleotide at said at least one attenuation locus in said PCR products with the corresponding nucleotide in the wild-type virus sequence, wherein the presence of the wild-type nucleotide at said at least one attenuation locus is indicative of the presence of a revertant virus in said sample.

2. The method of claim 1, wherein said attenuated virus is an RNA virus selected from dengue virus, poliovirus, rubella virus, measles virus, yellow fever virus, mumps virus and influenza virus.

3. The method of claim 1, wherein the at least one revertant virus is present in the sample containing the attenuated virus in a quantity of less than 1%, based on the total amount of virus in the sample.

4. The method of claim 1, wherein the at least one revertant virus differs from the attenuated virus by less than 10 nucleotides.

5. The method of claim 1, wherein said PCR products have a length of at least 3 kb.

6. The method of claim 1, wherein said PCR products comprise at least two attenuation loci.

7. The method of claim 1, wherein said at least one attenuation locus is selected from positions 57, 2579 and/or 5270 of the nucleotide sequence of dengue virus type 2 according to SEQ ID NO: 23.

8. The method of claim 1, wherein prior to step a) the live, attenuated virus is grown in Vero cells.

9. The method of claim 1, wherein the virus nucleic acids are isolated virus RNA and step b) comprises a step of reverse transcription of the isolated virus RNA in the presence of a reverse transcriptase using a reverse primer.

10. The method of claim 9, wherein the reverse transcription is performed using as reverse primer an oligonucleotide comprising the nucleotide sequence of SEQ ID NO:6 (3R) at a temperature in the range from about 45° C. to about 65° C.

11. The method of claim 1, wherein the PCR amplification reaction comprises the use of a first primer comprising the nucleotide sequence of SEQ ID NO:6 (3R) and a second primer comprising the nucleotide sequence of SEQ ID NO:9 (5F).

12. The method of claim 1, wherein prior to step d) the method additionally comprises:

- i) end repairing of the PCR products;
- ii) optionally, purifying the end repaired PCR product; and
- iii) connecting two complementary strands of the PCR product prepared in step c) with a linking oligonucleotide.

otide, wherein the 3' end of one complementary strand is linked to the 5' end of the other complementary strand via the linking oligonucleotide in each double stranded nucleic acid fragment, whereby the linking oligonucleotide provides a single stranded portion of the resulting linked nucleic acid fragments.

13. A method for the quantitative analysis of a mixture of virus haplotypes in a virus sample comprising the steps of:

- a) purifying virus nucleic acids from said sample;
- b) preparing double-stranded DNA templates from said virus nucleic acids by using a primer comprising a unique molecular identifier (UMI) sequence such that said double-stranded DNA templates comprise the UMI sequence;
- c) amplifying said double-stranded DNA templates by polymerase chain reaction (PCR) thereby generating PCR products having a length of at least 2 kb and comprising a unique molecular identifier (UMI) sequence;
- d) sequencing said PCR products by single-molecule real-time (SMRT) sequencing; and
- e) determining the relative amounts of said PCR products by analysing a UMI count and a Read count, thereby quantifying the virus haplotypes contained in the sample.

14. The method of claim **13**, wherein any of the virus haplotypes is present in the virus sample in a quantity of less than 1% based on the total amount of virus in the sample.

15. The method of claim **13**, wherein at least two of the virus haplotypes differ from each other by less than 10 nucleotides.

16. The method of claim **13**, wherein the virus sample comprises an RNA virus, selected from dengue virus, poliovirus, rubella virus, measles virus, yellow fever virus, mumps virus and influenza virus.

17. The method of claim **13**, wherein said virus sample comprises a live, attenuated dengue virus comprising at least one attenuation locus selected from positions 57, 2579 and/or 5270 of the nucleotide sequence of dengue virus type 2 according to SEQ ID NO: 23.

18. The method of claim **13**, wherein said PCR products have a length of at least 3 kb.

19. The method of claim **13**, wherein the virus nucleic acids are isolated virus RNA and step b) comprises a step of reverse transcription of the isolated virus RNA in the presence of a reverse transcriptase using a reverse primer additionally comprising at its 5' end a UMI sequence.

20. The method of claim **19**, wherein the reverse transcription is performed using as reverse primer an oligonucleotide comprising the nucleotide sequence of SEQ ID NO: 22 (UMI 3R) at a temperature in the range from about 45° C. to about 65° C.

21. The method of claim **13**, wherein the PCR amplification reaction comprises the use of a first primer comprising

the nucleotide sequence of SEQ ID NO: 21 (Primer IIA) and a second primer comprising the nucleotide sequence of SEQ ID NO: 9 (5F).

22. The method of claim **21**, wherein after the amplification step c) and prior to the sequencing step d) the method additionally comprises:

- i) end repairing of the PCR products;
- ii) optionally, purifying the end repaired PCR product; and
- iii) connecting the two complementary strands of the double stranded nucleic acid fragments prepared in step (b) with a linking oligonucleotide, wherein the 3' end of one complementary strand is linked to the 5' end of the other complementary strand via the linking oligonucleotide in each double stranded nucleic acid fragment, whereby the linking oligonucleotide provides a single stranded portion of the resulting linked nucleic acid fragments.

23. The method according to claim **1**, wherein said sample is a composition comprising an attenuated virus, or a sample from a patient.

24. The method of claim **1**, wherein prior to step d) the method comprises:

- i) determining the nucleotide sequence of the PCR products at at least one of the attenuation locus separately by Sanger sequencing;
- ii) comparing at least one nucleotide at said at least one attenuation locus with the wild-type virus sequence; and
- iii) if the wild-type sequence is determined at at least one attenuation locus, sequencing said PCR products by SMRT sequencing of step d) is carried out.

25. (canceled)

26. (canceled)

27. A kit for detecting the presence of at least one revertant virus in a sample containing an attenuated dengue virus comprising a first primer selected from an oligonucleotide having the nucleotide sequence of SEQ ID NO: 1, SEQ ID NO: 3, SEQ ID NO: 5, SEQ ID NO: 7, SEQ ID NO: 9 or SEQ ID NO: 11, and a second primer selected from an oligonucleotide having the nucleotide sequence of SEQ ID NO: 2, SEQ ID NO: 4, SEQ ID NO: 6, SEQ ID NO: 8, SEQ ID NO: 10 or SEQ ID NO: 12.

28. A kit for the quantitative analysis of a mixture of virus haplotypes in a dengue virus sample, comprising a first primer comprising the nucleotide sequence of SEQ ID NO: 9 (5F) and a second primer comprising the nucleotide sequence of SEQ ID NO: 21 (Primer IIA).

29. (canceled)

30. (canceled)

31. The kit of claim **27**, wherein the first primer comprises the nucleotide sequence of SEQ ID NO: 9 (5F) and the second primer comprises the nucleotide sequence of SEQ ID NO: 6 (3R).

* * * * *